Fast T Wave Detection Calibrated by Clinical Knowledge with Annotation of P and T Waves

Mohamed Elgendi^{1,2,*}, Bjoern Eskofier³ and Derek Abbott⁴

¹ Electrical and Computer Engineering in Medicine Group, University of British Columbia

 $^{\ 2}$ Department of Computing Science, University of Alberta

 3 Pattern Recognition Lab, Friedrich-Alexander University Erlangen-Nuernbeg

 4 School of Electrical and Electronic Engineering, University of Adelaide

July 21, 2015

Abstract

Background: There are limited studies on the automatic detection of T waves in arrhythmic electrocardiogram (ECG) signals. This is perhaps because there is no available arrhythmia dataset with annotated T waves. There is a growing need to develop numerically-efficient algorithms that can accommodate the new trend of battery-driven ECG devices. Moreover, there is also a need to analyze long-term recorded signals in a reliable and time-efficient manner, therefore improving the diagnostic ability of mobile devices and point-of-care technologies. Methods: Here, the T wave annotation of the well-known MIT-BIH arrhythmia database is discussed and provided. Moreover, a simple fast method for detecting T waves is introduced. A typical T wave detection method has been reduced to a basic approach consisting of two moving averages and dynamic thresholds. The dynamic thresholds were calibrated using four clinically known types of sinus node response to atrial premature depolarization (compensation, reset, interpolation, and reentry). Results: The determination of T wave peaks is performed and the proposed algorithm is evaluated on two well-known databases, the QT and MIT-BIH Arrhythmia databases. The detector obtained a sensitivity of 97.14% and a positive predictivity of 99.29% over the first lead of the validation databases (total of 221,186 beats). Conclusions: We present a simple yet very reliable T wave detection algorithm that can be potentially implemented on mobile battery-driven devices. In contrast to complex methods, it can be easily implemented in a digital filter design.

1 Introduction

According to the World Health Organization, cardiovascular diseases (CVDs) are the number one cause of death globally; more people die annually from CVDs than from any other cause [1]. An estimated 17.3 million people died from CVDs in 2008, representing 30% of all global deaths [1]. Of these deaths, an estimated 7.3 million were due to coronary heart disease and 6.2 million were due to stroke [1]. Thus, medical researchers have placed significant importance on cardiac health research. This has led to a strong focus on technological advances with respect to cardiac function assessment. One such research pathway is the improvement of conventional cardiovascular diagnosis technologies used in hospitals and clinics.

The most common clinical cardiac test is electrocardiogram (ECG) analysis as it is simple, risk-free, and inexpensive [2]. The signal of each heart beat contains five main waves: the P, Q, R, S, and T waves. The automatic detection of these waves is critical for reliable cardiovascular assessment, such as diagnosing cardiac arrhythmias [3, 4, 5, 6], understanding autonomic regulation of the cardiovascular system during sleep and hypertension [7, 8], detecting breathing disorders such as obstructive sleep apnea syndrome [9, 10], and monitoring other structural or functional cardiac disorders.

The detection of R peaks and QRS complexes has been extensively investigated over the past two decades [11]. Conversely, T wave detection has not been investigated as widely as QRS detection, and the T wave detection problem is still far from being solved [12, 13, 14, 15]. Reliable T wave detection is more challenging than QRS complex detection for several reasons, including low amplitudes, low signal-to-noise ratio (SNR), amplitude and morphology variability, and possible overlapping of the P wave and T wave [11, 16]. Nevertheless, accurate T wave detection is mandatory for a variety of (differential) diagnostic tasks, such as acute coronary syndrome [17], acute myocardial infarction [18], or potentially fatal arrhythmias [19]. To our knowledge, there was no attempt to develop an ECG detector based on existing clinical knowledge. Thus, in this paper, we investigated the possibility of applying clinical knowledge to build a T wave detection algorithm.

In addition, in the near future, it is expected that Holter devices, which are traditionally used for ECG analysis in the clinic, will be replaced by portable battery-operated devices, such as mobile phones [11]. Therefore, there is a need for a simple, fast, and computationally efficient algorithm to detect arrhythmias efficiently in real time. To develop fast robust algorithms for detecting Arrhythmia in ECG collected by portable, wearable, and battery-driven devices, first we require fully annotated arrhythmia ECG signals as a benchmark for evaluation. Unfortunately, the MIT-BIH Arrhythmia database [20] includes *only* the annotations of R peaks. Therefore, in this study, we annotated T waves in the MIT-BIH Arrhythmia database [21, 20]. Moreover, a new fast robust algorithm consisting of two moving averages that are calibrated by a clinical knowledge base is presented.

2 Materials and Methods 2.1 Data Used

Several standard ECG databases are available for the evaluation of QRS detection algorithms for ECG signals. Most of these databases contain annotated files for R peaks but not for T waves. To demonstrate the applicability of the algorithm presented in this paper, two databases are used in this study: one self-annotated database and one standard annotated database.

2.1.1 Database Annotated Under This Study

An expert manually annotated the P and T peaks of the MIT-BIH Arrhythmia database [21, 20] to be used in evaluation for the following reasons:

- 1. The MIT-BIH database contains 30-min recordings for each patient, which is considerably longer than the records in many other databases, such as the Common Standards for Electrocardiography database, which contains 10-s recordings [22].
- 2. The MIT-BIH Arrhythmia database contains records of normal ECG signals and records of ECG signals that are affected by non-stationary effects, low SNR, premature atrial complexes, premature ventricular complexes, left bundle blocks, and right bundle blocks. This provides an opportunity to test the robustness of T wave detection methods.
- 3. The database contains 23 records (the "100 series") that were chosen at random from a set of more than 4000 24-h Holter tapes, and 25 records (the "200 series") that were selected from the same set, including a variety of rare and clinically important ECG segments [20]. Several records in the 200 series have abnormal rhythms and QRS

morphologies and they suffer from a low SNR. These issues are expected to present significant difficulties for any ECG signal analysis algorithm [20].

Figures 1–3 demonstrate examples for annotation of T waves for different beats in the MIT-BIH Arrhythmia database. There was no automated aid provided during the annotation process and only channel one (Lead I) was annotated. For special cases, such as biphasic T waves, the middle point of the wave was considered as a T wave. The annotation file of P and T waves can be downloaded from [23].

2.1.2 Standard Annotated Database

As the MIT-BIH database is self-annotated, the validation of the detector must be carried out using a standard annotated database. For this purpose, the easily-available QT database [24] is used. This database was annotated by two cardiologists and includes different morphologies such as ST change, supraventricular arrhythmia, normal sinus rhythm, sudden death, and longterm ECG signals. The two cardiologists annotated only selected beats (3542 beats in a file called ".q1c") in all recordings except two recordings: "sel35" and "sel37". However, the automatic annotation of the whole database was carried using *ecgpuwave* software, which is saved in the ".pu" file. In this work, the T peaks of Lead I of the whole QT database are used for validation as they are more salient and certain compared to the onset and offset. In this work, the T peaks of the whole QT database are used for validation as they are more salient and certain compared to the onset and offset. Moreover, once the T peak is detected correctly, searching for the onset and offset is a relatively easy step; however, this is not the focus of this study.

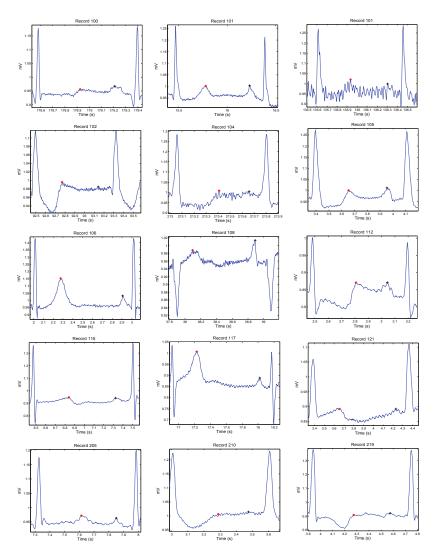


Figure 1: Annotation of P and T waves in normal beats. Here, "+" represents the P wave and "*" represents the T wave. Note that the higher peaks represent the QRS complex.

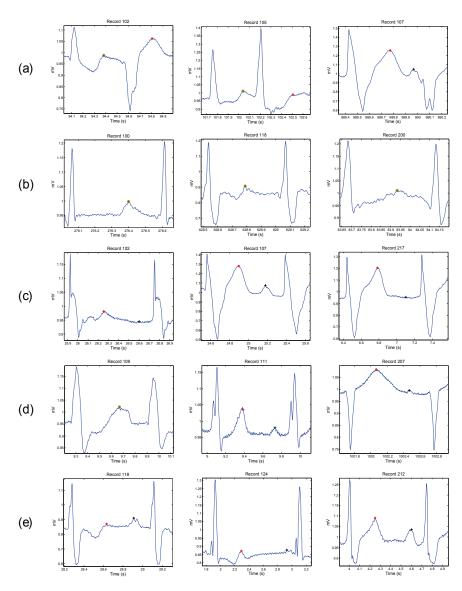


Figure 2: Annotation of Pand T waves in irregular heart beats. Each row contains three different morphologies for a certain type of arrhythmia: (a) premature ventricular beats; (b) premature atrial beats; (c) paced beats; (d) left bundle branch block beats; (e) right bundle branch block beats. Here, "+" represents the P wave and "*" represents the T wave, while the green circle with asterisk represents merged P and T waves.

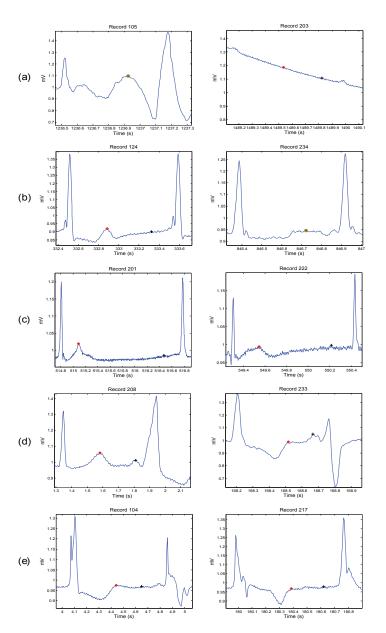


Figure 3: Annotation of P and T waves in unusual beats. Each row contains two different morphologies for a certain type of unusual beats: (a) unclassified beats; (b) nodal premature beat; (c) nodal escape beat; (d) fusion of ventricular and normal beat; (e) fusion of paced and normal beat. Here, "+" represents the P wave and "*" represents the T wave, while the green circle with asterisk represents merged P and T waves.

2.2 T Wave Detection Algorithm

In this study, a fast robust knowledge-based T wave detection algorithm is discussed and evaluated. The algorithm is based on the framework proposed by Elgendi for detecting QRS complexes in ECG signals [25, 26], for detecting systolic waves in photoplethysmogram signals [27], detecting a waves [28], and detecting c, d, e waves [29] in acceleration photoplethysmogram (PPG) signals. We build upon this approach to detect T waves. The method consists of three main stages: pre-processing (clinical knowledge, bandpass filtering, squaring, and QRS removal), feature extraction (generating potential blocks using two moving averages), and classification (thresholding). The structure of the algorithm is given in Figure 4.

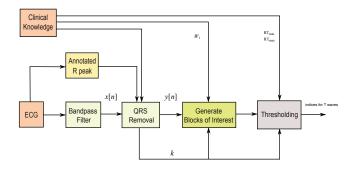


Figure 4: Structure of the T waves detection algorithm. The algorithm consists of three main stages: pre-processing (clinical knowledge, bandpass filtering, squaring, and QRS removal), feature extraction (generating potential blocks using two moving averages), and classification (thresholding).

2.2.1 Bandpass Filter

Most of the energy of T waves lies below 10 Hz [30, 31]; thus, a zero-phase second-order Butterworth filter, with a 0.5–10 Hz bandpass, is implemented to remove the baseline wander and high frequencies that do not contribute to the T waves. The output of the zero-phase Butterworth filter applied to

the ECG signal produces a filtered signal x[n].

2.2.2 QRS Removal

The QRS removal step based on the relative distance of the T wave to its associated R peak. Removing the QRS complex has two advantages: (1) T waves become the dominant feature in the processed signal; and (2) it simplifies the search for T waves relative to the position of R peaks. In this study, as proof of concept, the R peaks provided in the MIT-BIH Arrhythmia and QT databases are used. Removing the QRS complex is performed by setting the signal to zero for the duration of the QRS complex.

The limits (thresholds) of the RT distance are determined using ECG clinical knowledge. The signal y[n] is initialized as equal to the filtered ECG x[n] signal. The QRS removal length thresholds are determined based on the clinical phases of the RR intervals. As the duration of the QRS complex varies with the heart beat type, a clinical database is required to remove the QRS complex, according to its type. Roskamm and Csapo divided the ECG into four phases: compensation, reset, interpolation, and reentry [32]. Based on their analysis, in the compensation phase (Figure 5a), the second beat (850 ms) is followed by a prolonged beat (1150 ms) to compensate the two beats duration (2000 ms). During the reset phase (Figure 5b), the second beat (650 ms) is followed by a prolonged beat (1150 ms), while in the interpolation (Figure 5c), the second beat (400 ms) is followed by an irregular beat (600 ms). Finally, in the reentry phase (Figure 5d), the second beat (300 ms) is followed by a rapid irregular beat (400 ms); however, an extra category is added to capture complex arrhythmias (repetitive, bigeminy, or trigeminy). The output of this stage will produce an updated y[n] signal.

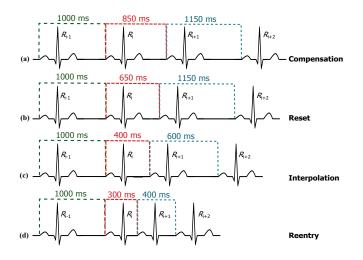


Figure 5: Types of sinus node response to atrial premature depolarization adapted from [32].

Based on the clinical information presented in Figure 5, the ratio of the second beat duration to the first and third beat duration generates a rule-based knowledge representation, as shown in Figure 6. For example, in the case of reentry (left branch of the flowchart in Figure 6), the first threshold is calculated by dividing the duration between the third and the fourth beat $(R_{i+2} - R_{i+1} = 400 \text{ ms})$ by the duration between the two and the third beat $(R_{i+1} - R_i = 300 \text{ ms})$ resulting as 1.33; however, this value was calculated with reference to the duration between first beat and second beat $(R_i - R_{i-1} = 1000 \text{ ms})$ —this generates the first rule: RR2 \leq 1.33 RR1, where $RR1 = R_i - R_{i-1}$, $RR2 = R_{i+1} - R_i$, and *i* is the beat index. The second threshold is the expected total duration of the second and third beats (300 ms + 400 ms), which equals 0.7 s—this generates the second rule: $(RR1 + RR2) \leq 0.7 f_s$. On the right side the flowchart in Figure 6, in the case of reentry, the first threshold is calculated by dividing the duration between the third and the fourth beat $(R_{i+1} - R_i = 300 \text{ ms})$ by the duration between the first and the second beat $(R_i - R_{i-1} = 1000 \text{ ms})$

resulting as 0.3—this generates the first rule: RR2 \leq 0.3 RR1, while the second threshold is based on the duration between the first and the second beat $(R_i - R_{i-1} = 1000 \text{ ms})$ which is 1 s (or f_s)—this generates the second rule: RR1 \leq f_s . Similarly, all other rules were created for the other four ECG phases.

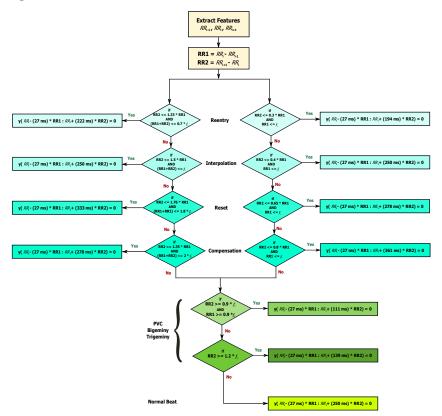


Figure 6: Rule-based knowledge representation of QRS removal based on the clinical knowledge shown in Figure 5.

During the QRS removal, the RR interval that satisfied each category is saved and referred to as RR_k , where k is the category type (compensation, reset, interpolation, reentry, and complex arrhythmias). The normalized RR intervals average in each category is calculated as $M_c = (\sum_{j=1}^{l} \operatorname{RR}_{c,j})/(lf_s)$, where l is the number of RR intervals saved in category c, and f_s is the sampling frequency. Dividing M_c by f_s is sufficient as it is equivalent to the average RR interval in healthy subjects.

2.2.3 Generating Blocks of Interest

Blocks of interest are generated using two event-related moving averages that demarcate the areas of T waves, a method which was first introduced in [33]. The particular method used to generate blocks of interest has been mathematically shown to detect a waves [28], QRS complexes [25], and systolic waves in PPG signals [27]. In this procedure, the first moving average (MA_{peak}) is used to emphasize the peak of the T wave area, as the dotted signal shown in Figure 7, and is given by

$$MA_{peak}[n] = \frac{1}{W_1} (y[n - (W_1 - 1)/2] + \dots + y[n] + \dots + y[n + (W_1 - 1)/2])$$
(1)

where W_1 represents the window size of approximately the peak duration of the T wave in ECG signals. The initial value for W_1 of 70 ms is determined by Trahanias [34]. However, as the ECG signals may contain different arrhythmias the value of W_1 will be calculated relative to the most frequent RR intervals in all five categories ($k = \max_c M_c$). Then, the value of $W_1 = (70 \text{ ms} \times f_s) \times k$, and the result is rounded to the nearest odd integer. The second moving average (MA_{Twave}) is used to emphasize the T wave area to be used as a threshold for the first moving average, shown as a dashed signal Figure 7, and is given by

$$MA_{Twave}[n] = \frac{1}{W_2} (y[n - (W_2 - 1)/2] + \dots + y[n] + \dots + y[n + (W_2 - 1)/2])$$
(2)

where W_2 represents the window size of approximately the T wave duration. The initial value for W_2 of 140 ms is determined by Laguna *et al.* [35]. However, as the ECG signals may contain different arrhythmias, the value of W_2 will be calculated relative to the most frequent RR intervals in all five categories (k). Then, the value of $W_2 = (140 \text{ ms} \times f_s) \times k$, and the result is rounded to the nearest odd integer. For example, the total values of W_1 and W_2 for detecting T waves in record 100 from MIT-BIH Arrhythmia database were 20 samples (55.6 ms) and 40 samples (111.2 ms); respectively.

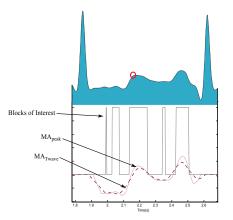


Figure 7: Demonstrating the effectiveness of using two moving averages to detect T waves. The dotted line is the first moving average, while the dashed line is the second moving average. The red circle is the detected T wave peak.

2.2.4 Thresholding

In this stage, the blocks of interest are generated by comparing the MA_{peak} signal with MA_{Twave} . Many blocks of interest will be generated, some of which will contain the T wave and others will contain P waves, U waves, and noise. Therefore, the next step is to reject blocks that result from noise. Rejection is based on the relative positions of P and T waves to R peaks and anticipated peak width.

To determine whether the detected blocks contain T waves or not, the number of blocks in each consecutive RR interval is counted. A threshold based on the distance of the maximum point within a block to the R peak is applied to distinguish P waves from T waves and noise, as shown in Figure 8. The search regions for T waves in terms of time occurrence with respect to the current R peak (R_i) and the next R peak (R_{i+1}) are calculated as

$$R_i T_{\min} = D_{\min} R_i R_{i+1} \tag{3}$$

$$\mathbf{R}_i \mathbf{T}_{\max} = \mathbf{D}_{\max} \mathbf{R}_i \mathbf{R}_{i+1} \tag{4}$$

where $R_i T_{min}$ represents the minimum dynamic interval between the T wave and the current R peak, $R_i T_{max}$ represents the maximum dynamic interval between the T wave and the current R peak, while $R_i R_{i+1}$ represents the interval between R_i and R_{i+1} . The exact values for D_{min} and D_{max} are 170 ms and 800 ms, respectively, as determined by Schimpf *et al.* [36] to represent the minimum RT durations for subjects with arrhythmia and maximum RT duration for healthy subjects. All detected blocks go through a durational threshold to reject the undesired blocks called THR₁, which rejects the blocks that contain P wave, U wave, and noise. By applying the THR₁ threshold, the accepted blocks will contain T peaks only,

$$THR_1 = W_1 \tag{5}$$

After applying the relative-position thresholds, there are three possibilities for the number of detected blocks within the area of interest:

- Zero: if there is no block detected, it means the algorithm failed to detect a T wave in the current RR interval.
- 2. **One:** if there is one detected block, it means the algorithm succeeds in detecting T wave, P and T waves are most likely merged within one block, which is marked as a circle with a black asterisk inside (see Figure 9i,j).

3. More than one: if there are multiple detected blocks, as shown in Figure 7, it means one of the detected blocks contains T waves. However, in this work the nearest block to the current R peak is considered a T wave.

The last stage is to find the maximum absolute value within each block to detect the peak of T wave. The detected T wave peaks are compared to the annotated T wave peaks to determine whether they were detected correctly. The search range for the true T wave peak is fixed to ± 30 ms for both databases, to ensure consistency of comparison. The search region of 30 ms is good enough for diagnostics as it is less than W_1 .

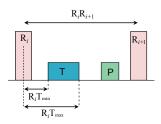


Figure 8: Search regions for T waves in terms of time occurrence with respect to the current R peak (R_i) and the next R peak (R_{i+1}). Where $R_i T_{\min}$ represents the minimum interval between the T wave and current R peak and $R_i T_{\max}$ represents the maximum interval between the T wave and the current R peak.

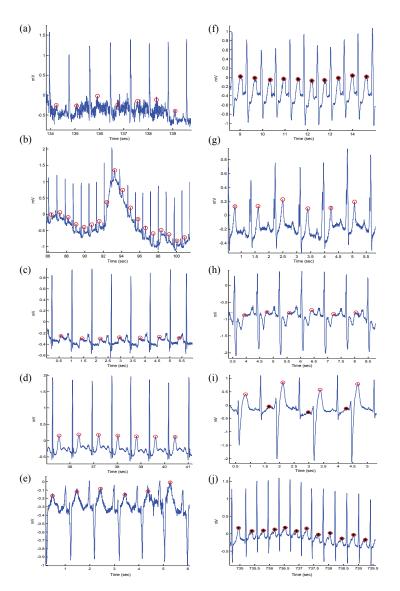


Figure 9: Demonstrating the performance of the proposed T wave detection algorithm on the MIT-BIH Arrhythmia database. The algorithm succeeds to detect T wave peaks in electrocardiogram (ECG) signals that contain: (a) high-frequency noise; (b) baseline wander; (c) normal sinus rhythm without U waves; (d) normal sinus rhythm with U waves; (e) normal sinus rhythm with negative polarization; (f) left bundle branch block (LBBB) beats with merged P and T waves; (g) LBBB beats; (h) right bundle branch block (RBBB) beats from record 118; (i) premature ventricular contraction (PVC) beats from record 200; (j) premature atrial contraction (PAC) beats from record 209. Here, the empty red circle represents the detected T wave while the a circle with a black asterisk represents detection of merged P and T waves.

3 Results

The algorithm was evaluated using the MIT-BIH database. The T waves were detected successfully even when the P and T waves are merged in Arrhythmia ECG signals that are affected by: high-frequency noise, baseline wander, normal sinus rhythm (NSR), left bundle branch block (LBBB), right bundle branch block (RBBB), premature ventricular contraction (PVC), and premature atrial contraction (PAC). All of the reasons for detection failure are described below. High-frequency noise results from the instrumentation amplifiers, recording system, and ambient electromagnetic signals received by the cables. The signal shown in Figure 9a has been corrupted by powerline interference at 60 Hz and its harmonics and other high frequencies. It can be seen that the proposed algorithm is robust to noise. Moreover, the proposed algorithm is not sensitive to baseline wander and detected the T waves correctly, as shown in Figure 9b. This is because the moving averages were applied to the bandpass filtered ECG signal—which is discussed in the Bandpass Filter Subsection.

The NSR is a normal ECG cycle; it is initiated by the sinoatrial node and consists of a P wave followed, after a brief pause, by a QRS complex and then a T wave [37]. The proposed algorithm correctly detected T waves in three types of normal beats: (1) NSR without U waves (record 100 of the MIT-BIH database), as shown in Figure 9c; (2) NSR with U waves (record 103), as shown Figure 9d; and (3) NSR with negative polarization (record 108), as shown Figure 9e. The LBBB results from conduction delays or blocks at any site in the intraventricular conduction system, including the main LBBB and the bundle of His. The result of an LBBB is an extensive reorganization of the activation pattern of the left ventricles [37]. The proposed algorithms successfully detected normal and merged P and T waves in two types of LBBBs: (1) LBBB beats with merged P and T waves (record 109), as shown in Figure 9f and (2) LBBB beats with normal T waves (record 111), as shown in Figure 9g. However, RBBB is a result of a conduction delay in a portion of the right-sided intra-ventricular conduction system. The delay can occur in the main RBBB itself, in the bundle of His, or in the distal right ventricular conduction system. The RBBBs may be caused by a minor trauma, such as right ventricular catheterization [37]. As shown in Figure 9h, the proposed algorithms succeeded in detecting the T waves in ECG signals of RBBB (record 118).

The PVCs are characterized by the premature occurrence of a QRS complex that is abnormal in shape, and that has a longer duration than normal QRS complexes, generally exceeding 120 ms [37]. The T wave is commonly large and opposite in direction to the major deflection of the QRS. In general, the QRS complex is not preceded by a P wave, but it can be preceded by a non-conducted sinus P wave occurring at the expected time [37]. In Figure 9i, a special case of PVC is shown, called bigeminy, where the premature ventricular beats occur after every normal beat in an alternating pattern.

The proposed algorithm succeeded in detecting the T waves in the normal beats and the T waves in the premature ventricular beats (record 200). Note that PACs are among the most common causes of irregular pulses and can originate from any area in the heart [37]. The impulse is discharged prematurely by an irritable focus in the atria giving rise to a distorted P wave, usually superimposed on the preceding T wave. As shown in Figure 9j, the proposed algorithms detected the merged T waves in PACs (record 209). As illustrated in Figure 9, the proposed method successfully detected T waves in ECG signals with a low SNR, baseline wander, and various arrhythmias. The performance of the T wave detection algorithms is evaluated using two statistical measures: SE = TP/(TP + FN) and +P = TP/(TP + FP), where TP is the number of true positives (T wave peak detected within the range of 30 ms of the annotated T wave peak), FN is the number of false negatives (annotated T wave peak has not been detected), and FP is the number of false positives (T wave peak detected outside the range of 30 ms of the annotated T wave peak). The sensitivity SE reports the percentage of true beats that were correctly detected by the algorithm. The positive predictivity +Preports the percentage of beat detections that were true beats.

The abnormal heart rhythms caused a large number of FNs compared to the FPs. Table 1 shows the result of T wave detection over 48 records of the MIT-BIH database. FNs are mainly caused by noise and PVC, as in record 219, and atrial fibrillation, as in record 202. The algorithm achieved a sensitivity of 99.86% and a positive predictivity of 99.65%, which are promising results for handling the non-stationary effects, low SNR, PACs, PVCs, LBBBs, and RBBBs in ECG signals.

3.1 Comparison of Performance on the QT Database

The detection performance on the QT database obtained by the proposed T wave detector record by record performance is shown in Tables 2 and 3. The overall comparison of our results with the existing T wave detection algorithms on the QT database is demonstrated in Table 4. It summarizes the performances in terms of number of beats, methodology, SE, and +P. Note that the proposed algorithm scored slightly higher overall performances (average of SE and +P) than Martínez *et al.* [38] and Laguna *et al.* [39] over the manually annotated T waves. It is clear that the proposed algorithm scored slightly higher overall performance over the manual of the performance over the performance performance over the performance over the performance performance over the performance over the performance performance over the performance p

111,201 automatically annotated heart beats. Moreover, the proposed T wave detector has not been re-tuned over any databases, thus the results are promising, and the algorithm can detect T peaks over different databases, sampling frequencies, types of arrhythmias, and noise.

3.2 Processing Time

Less computational time is achieved when the simplest method is used (*i.e.*, the algorithm presented in this paper requires less computational time). This is advantageous in terms of future development of wearable and portable diagnostic devices, and in terms of helping online and real-time diagnoses.

It is misleading to suggest that mentioning the average speed of the proposed detector, over a certain time length of ECG signal, would provide a comparative result. This is because the processing time depends on the number of beats within each ECG recording, not on the record length. In this study, the T wave detector was implemented in MATLAB 2010b (The MathWorks, Inc., Natick, MA, USA) on IntelTM i5 CPU 2.27 GHz.

It is worth noting that the number of beats of the 15-min recordings category in QT database was relatively consistent—with a mean \pm SD, number of beats 1059 \pm 275—over all records of this category. On the contrary, the 30-min beat average, in MIT-BIH database, was 2291 with an SD of 448 beats. As the processing time depends on the number of beats rather than the recording length [25], we found, for example, the proposed detector took 2.58 s to process record 215-MITDB contains 3362 beats, while it took 0.28 s to process 527 beats in record 33-QTDB. In general, without taking the number of beats into consideration, the speed of the proposed detector is fast. The suggested detector handles 15-min recordings in 0.52 \pm 0.14 s, while it takes 1.98 \pm 0.32 s to handle 30-min ECG recordings.

Table 1: T wave peak detection performance over the annotated MIT-BIH Arrhythmia database [21, 20]. To evaluate the performance of the T wave detection algorithm, two statistical measures are used: SE = TP/(TP + FN) and +P = TP/(TP + FP), where TP is the number of true positives (T wave peak detected within the range of 30 ms of the annotated T wave peak), FN is the number of false negatives (annotated T wave peak has not been detected), and FP is the number of false positives (T wave peak detected outside the range of 30 ms of the annotated T wave peak detected outside the range of 30 ms of the annotated T wave peak.

Record	No of beats	TP	FP	FN	SE (%)	+P (%)
100	2274	2272	0	0	100.00	100.00
101	1866	1863	1	4	99.79	99.95
102	2187	2185	0	0	100.00	100.00
103	2084	2082	0	4	99.81	100.00
104	2229	2227	0	1	99.96	100.00
105	2602	2586	0	2	99.92	100.00
106	2026	2024	0	56	97.23	100.00
107	2136	2134	0	3	99.86	100.00
108	1763	1757	0	13	99.26	100.00
109	2533	2530	0	0	100.00	100.00
111	2123	2121	0	16	99.25	100.00
112	2539	2537	0	0	100.00	100.00
113	1794	1792	0	0	100.00	100.00
114	1890	1885	0	69	96.34	100.00
115	1953	1951	0	19	99.03	100.00
116	2395	2392	0	2	99.92	100.00
117	1535	1533	0	0	100.00	100.00
118	2278	2276	0	4	99.82	100.00
119	1988	1986	0	4	99.80	100.00
121	1863	1860	0	46	97.53	100.00
122	2476	2474	0	0	100.00	100.00
123	1519	1517	0	0	100.00	100.00
124	1619	1617	0	7	99.57	100.00
200	2601	2599	0	9	99.65	100.00
201	1949	1947	0	57	97.07	100.00
202	2138	2134	0	113	94.70	100.00
202	2988	2965	0	1	99.97	100.00
205	2656	2556	0	0	100.00	100.00
205	2324	2139	0	9	99.58	100.00
207	2953	2139 2949	0	9	100.00	100.00
208	3006	3003	0	5	99.83	100.00
210	2652	2637	0	0	99.85 100.00	
210						100.00
212	2748	2746	0	0	100.00	100.00
	3250	3247	0	0	100.00	100.00
214 215	2262	2184	0	0	100.00	100.00
	3362	3354	0	0	100.00	100.00
217	2208	2205	0	3	99.86	100.00
219	2154	2152	0	144	93.31	100.00
220	2048	2046	0	2	99.90	100.00
221	2427	2424	0	0	100.00	100.00
222	2485	2472	0	33	98.67	100.00
223	2604	2601	0	1	99.96	100.00
228	2060	2056	0	52	97.47	100.00
230	2256	2254	0	39	98.27	100.00
231	1571	1569	0	0	100.00	100.00
232	1783	1781	0	1	99.94	100.00
233	3077	2914	0	1	99.97	100.00
234	2751	2749	0	0	100.00	100.00

Table 2: T wave peak detection performance over the manually annotated 11 recordings of the QT database [24]. To evaluate the performance of the T wave detection algorithm, two statistical measures are used: SE = TP/(TP + FN) and +P = TP/(TP + FP), where TP is the number of true positives (T wave peak detected within the range of 30 ms of the annotated T wave peak), FN is the number of false negatives (annotated T wave peak has not been detected), and FP is the number of false positives (T wave peak detected outside the range of 30 ms of the annotated T wave peak).

Record	No of beats	TP	FP	FN	SE (%)	+P (%)
sel100	30	30	0	0	100.00	100.00
sel102 sel103	85 30	85 30	0	0	100.00 100.00	100.00 100.00
sel103	77	75	0	0	100.00	100.00
sell14	50	50	0	0	100.00	100.00
sel116	50	49 30	0	0	100.00	100.00
sel117 sel123	30 30	30	0	0	100.00	100.00
sel213	71	71	0	0	100.00	100.00
sel221	30	29	1	1	96.67	96.67
sel223	31	31	0	0	100.00	100.00
sel230 sel231	50 50	42 50	9	8 0	84.00 100.00	82.35 100.00
sel232	30	30	0	0	100.00	100.00
sel233	30	30	0	0	100.00	100.00
sel301	30	30	0	0	100.00	100.00
sel302 sel306	30 36	30 32	0	0	100.00	100.00 100.00
sel300	30	30	0	0	100.00	100.00
sel308	50	40	10	10	80.00	80.00
sel310	30	30	0	0	100.00	100.00
sel803 sel808	30 30	30 30	4	0	100.00	88.24 100.00
sel808 sel811	30	30	0	0	100.00	100.00
sel820	30	30	0	0	100.00	100.00
sel821	30	30	ō	0	100.00	100.00
scl840	70	70	0	0	100.00	100.00
sel847	33	33	0	0	100.00	100.00
sel853 sel871	30 70	30 70	0	0	100.00	100.00 100.00
sel871 sel872	70 30	70	0	0	100.00	100.00
sel873	33	33	0	0	100.00	100.00
sel883	30	30	0	0	100.00	100.00
sel891	71 30	71 30	0	0	100.00	100.00 100.00
sel16265 sel16272	30 30	30 30	0	0	100.00	100.00
sel16272 sel16273	30	30	0	0	100.00	100.00
sel16420	30	30	0	0	100.00	100.00
sel16483	30	30	0	0	100.00	100.00
sel16539	30	30	0	0	100.00	100.00
sel16773 sel16786	30 30	21 30	9	9	70.00 100.00	70.00 100.00
sel16795	30	30	0	0	100.00	100.00
sel17453	30	30	0	0	100.00	100.00
sele0104	30	30	0	0	100.00	100.00
sele0106	30	30	0	0	100.00	100.00
sele0107 sele0110	34 30	34 30	0	0	100.00	100.00 100.00
selc0110	30	30	0	0	100.00	100.00
sele0112	50	50	0	0	100.00	100.00
sele0114	30	30	ō	ō	100.00	100.00
sele0116	30	30	0	0	100.00	100.00
sele0121 sele0122	30 30	30 30	0	0	100.00	100.00 100.00
sele0122 sele0124	30	50	0	0	100.00	100.00
sele0126	30	25	4	5	83.33	86.21
sele0129	30	30	0	0	100.00	100.00
sele0133	30	30	0	0	100.00	100.00
sele0136 sele0166	30 36	30 36	0	0	100.00 100.00	100.00 100.00
sele0100	30	30	0	ő	100.00	100.00
sele0203	30	30	ō	0	100.00	100.00
sele0210	30	30	0	0	100.00	100.00
sele0211	30	30	ō	ō	100.00	100.00
sele0303 sele0405	30 30	30 30	0	0	100.00	100.00 100.00
sele0405 sele0406	30	30	0	0	100.00	100.00
sele0409	30	30	ō	ò	100.00	100.00
sele0411	30	30	0	0	100.00	100.00
sele0509 sele0603	30 30	30 30	0	0	100.00	100.00
sele0603 sele0604	30 30	30 30	0	0	100.00	100.00
sele0606	30	30	ò	ò	100.00	100.00
sele0607	30	30	0	0	100.00	100.00
sele0609	30	30	0	0	100.00	100.00
sele0612 sele0704	30 30	30 30	0	0	100.00 100.00	100.00 100.00
sele0704 sel30	30 30	30 30	0	0	100.00	100.00 100.00
sel30 sel31	30	30	0	4	100.00 86.67	100.00 86.67
sel32	30	30	0	0	100.00	100.00
sel33	30	30	ō	ō	100.00	100.00
sel34	30	30	0	0	100.00	100.00
sel36 sel38	31 30	30 30	1	1	96.77	96.77 100.00
sel38 sel40	30	30	0	0	100.00	100.00
sel41	30	30	0	ò	100.00	100.00
sel42	30	30	0	0	100.00	100.00
sel43	30	28	0	0	100.00	100.00
sel44	30	30	0		100.00	100.00
sel45 sel46	30 30	29 30	1	1	96.67 100.00	96.67 100.00
sel46 sel47	30 30	30 27	0	0	100.00 93.10	100.00 93.10
sel48	30	30	0	0 0	100.00	100.00
sel49	30	30	0	0	100.00	100.00
sel50	30	30	ō	ō	100.00	100.00
sel51	32 30	32 30	0	0	100.00	100.00
sel51 sel52	30 30	30 30	0	0	100.00	100.00
sel52 sel17152	30	30	0	0	100.00	100.00
sel14046	31	31	0	0	100.00	100.00
		30	ō	0	100.00	100.00
sel14157 sel14172	30 50	30	0	0	100.00	100.00

Table 3: T wave peak detection performance over the automatically annotated QT database [24]. To evaluate the performance of the T wave detection algorithm, two statistical measures are used: SE = TP/(TP + FN) and +P = TP/(TP + FP), where TP is the number of true positives (T wave peak detected within the range of 30 ms of the annotated T wave peak), FN is the number of false negatives (annotated T wave peak has not been detected), and FP is the number of false positives (T wave peak detected outside the range of 30 ms of the annotated T wave peak).

beaks beaks atlog 1134 atlog 1088 atlog 1081 atlog 1081 atlog 1081 atlog 1031 atlog 1032 atlog 1032 atlog 1032 atlog 1032 atlog 1032 atlog 1032 atlog 1031 atlog 1031 atlog 1032 atlog 1031 atlog 1032 atlog 1031 atlog 1031 atlog 1031 atlog 1031 atlog 1032 <	1132 1136 1036 1046 1047 1141 1541 1541 1541 1541 1541 1547 1541 1547 1541 1547 1567 1567 1573 1568 1551 1552 1553 1517 1553 1517 1553 1115 1115 1115 1115 1115 1115 1115 1115 1115 1115 1115 1116 11178 1110 11110 11110 11111 11110 11110 11110 11110 11110 11110	$ \begin{smallmatrix} 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0$	1 2 5 0 7 25 1 0 2 1 0 2 1 0 2 1 1 2 1 1 2 30 1 1 2 3 2 3 3 2 1 1 2 3 3 2 1 1 2 1 1 1 2 3 3 2 1 1 1 2 3 2 3 2	9091 9082 9082 9082 9082 9082 9082 9082 9082	100.00 100.00 99.63 99.85 100.00 99.63 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 98.97 100.00 98.93 100.00 98.93 100.00 98.72 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.46 100.00 99.47 100.00 99.46 100.00 99.46 100.00
eld03 1048 eld04 1047 eld04 1107 eld14 1107 eld14 1107 eld14 1107 eld14 1107 eld14 1107 eld13 1107 eld14 1107 eld12 1101 eld21 1111 eld21 1111 eld21 1111 eld21 1111 eld21 1111 eld21 1111 eld213 1111 eld214 1111 eld215 1111 eld216 1111 eld217 1111 eld218 1111 eld211 11111 eld211	1046 1107 114 114 114 114 114 114 114 114 114 11	$\begin{smallmatrix} 4 & 9 & 1 \\ 0 & 0 & 0 & 0 \\ 0 & 0 & 0 & 0 \\ 0 & 0 &$	5 10 7 7 25 1 0 2 116 6 2000 1 19 1 20 1 2 1 0 2 3 0 1 2 1 3 84 9 1 2 1 3 8 4 9 1 2 1 3 8 4 1 1 1 5 1 5 1 5 1 5 1 5 1 5 1 5 1 5 1 5	$\begin{array}{c} 0.9, 22, \\ 0.9, 22, \\ 0.9, 10 \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 87, \\ 0.9, 10, \\ 0.9, 10, \\ 0.9, 10, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, 11, \\ 0.9, 10, \\ 0.9, $	99.62 99.18 99.19 99.00.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 99.70 100.00 99.53 100.00 99.53 100.00 99.53 100.00 99.53 100.00 99.55 100.000
elibei 1109 elibei 1101 elibei 1111 elibei </td <td>1107 864 164 164 164 164 164 164 164 164 164 1</td> <td>9 1 0 0 0 1 0 0 1 0 0 0 0 0 0 0 0 0 0 0</td> <td>10 7 215 1 1 6 2 10 6 2 00 1 12 2 3 8 4 2 3 8 4 2 3 2 2 3 3 2 2 1 5 1 1 1 2 1 5 1 1 2 1 6 6 00 1 1 1 2 2 3 3 2 2 3 8 3 2 2 1 1 6 6 0 0 2 1 1 6 6 0 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 2 3 3 2 2 3 2 3 2 2 3 2 2 1 1 6 6 0 0 1 1 1 2 2 3 3 2 2 3 2 2 1 1 6 6 0 0 1 1 1 2 2 3 3 2 2 1 1 1 1 2 2 3 3 2 2 1 1 1 2 1 1 2 1 1 1 2 2 3 3 2 2 1 1 1 2 2 3 3 2 2 1 1 1 1</td> <td>$\begin{array}{c} 99, 10\\ 99, 10\\ 99, 10\\ 90, 10\\$</td> <td>$\begin{array}{c} 99.19\\ 99.88\\ 100.60\\ 1$</td>	1107 864 164 164 164 164 164 164 164 164 164 1	9 1 0 0 0 1 0 0 1 0 0 0 0 0 0 0 0 0 0 0	10 7 215 1 1 6 2 10 6 2 00 1 12 2 3 8 4 2 3 8 4 2 3 2 2 3 3 2 2 1 5 1 1 1 2 1 5 1 1 2 1 6 6 00 1 1 1 2 2 3 3 2 2 3 8 3 2 2 1 1 6 6 0 0 2 1 1 6 6 0 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 1 5 6 0 0 1 1 1 2 2 3 3 2 2 3 2 3 2 2 3 2 2 1 1 6 6 0 0 1 1 1 2 2 3 3 2 2 3 2 2 1 1 6 6 0 0 1 1 1 2 2 3 3 2 2 1 1 1 1 2 2 3 3 2 2 1 1 1 2 1 1 2 1 1 1 2 2 3 3 2 2 1 1 1 2 2 3 3 2 2 1 1 1 1	$\begin{array}{c} 99, 10\\ 99, 10\\ 99, 10\\ 90, 10\\$	$\begin{array}{c} 99.19\\ 99.88\\ 100.60\\ 1$
atti 4 S67 atti 4 S67 atti 5 S62 atti 5 S62 </td <td>864 184 763 763 124 1307 124 1307 864 124 1307 864 125 126 87 126 87 126 885 1208 821 208 901 902 903 904 905</td> <td>$1 \\ 0 \\ 0 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 1 \\$</td> <td>25 0 2 116 6 2 2 19 19 112 0 2 2 0 1 112 0 2 2 3 1 12 12 3 8 4 2 1 3 3 2 2 3 3 2 2 1 1 1 12 0 2 3 0 1 12 12 0 2 3 12 13 8 42 12 12 0 2 13 14 5 12 12 0 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 3 2 12 12 12 2 3 2 3 2 2 3 2 3 2 2 1 12 12 12 12 13 13 13 14 12 12 12 12 12 12 13 13 14 14 12 12 13 13 14 14 12 12 13 13 14 14 14 12 14 14 14 12 14 14 14 14 14 14 14 14 14 14 14 14 14</td> <td>$\begin{array}{c} 00, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10$</td> <td>99.88 100.00 1100.00 99.94 100.00 99.94 100.00 99.94 100.00 97.92 97.92 97.92 100.00 97.92 100.00 97.93 100.00 99.93 100.00 99.93 100.00 99.93 100.00 99.93 100.00 99.95 100.00 100.00 99.95 100.00 100.</td>	864 184 763 763 124 1307 124 1307 864 124 1307 864 125 126 87 126 87 126 885 1208 821 208 901 902 903 904 905	$1 \\ 0 \\ 0 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 1 \\ $	25 0 2 116 6 2 2 19 19 112 0 2 2 0 1 112 0 2 2 3 1 12 12 3 8 4 2 1 3 3 2 2 3 3 2 2 1 1 1 12 0 2 3 0 1 12 12 0 2 3 12 13 8 42 12 12 0 2 13 14 5 12 12 0 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 12 12 0 2 2 3 2 12 12 12 2 3 2 3 2 2 3 2 3 2 2 1 12 12 12 12 13 13 13 14 12 12 12 12 12 12 13 13 14 14 12 12 13 13 14 14 12 12 13 13 14 14 14 12 14 14 14 12 14 14 14 14 14 14 14 14 14 14 14 14 14	$\begin{array}{c} 00, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10, 10\\ 0, 10$	99.88 100.00 1100.00 99.94 100.00 99.94 100.00 99.94 100.00 97.92 97.92 97.92 100.00 97.92 100.00 97.93 100.00 99.93 100.00 99.93 100.00 99.93 100.00 99.93 100.00 99.95 100.00 100.00 99.95 100.00 100.
eli17 766 1131 754 1131 754 1131 754 1131 754 1131 754 1131 754 1131 754 1131 754 1131 754 1232 1309 1232 1309 1232 1309 1323 1552 1324 1324 1324 1324 1324 1324 1324 1324 1324 1324 1324 1324 1334 1532 1340 1324 1340 1324 1340 1437 1340 1311 1340 1311 1340 1311 1340 1311 1340 1311 1340 1311 1340 1311 1340 1311 1340	764 714 124 1307 721 1307 723 725 725 725 725 725 725 725 725 725 725	0 0 1 1 0 1 1 0 1 1 0 1 1 0 0 1 1 0 0 1 1 0 0 1 1 0 0 1 1 0 0 1 1 0 0 1 1 0 0 0 1 1 0 0 0 1 1 0 0 0 1 1 0 0 0 1 0 0 0 1 0	1 0 2 116 6 200 119 112 0 2 2 30 1 2 1 3 3 3 3 3 3 3 3 3 3 3 3 3 3 3 3 3	$\begin{array}{c} 99.87\\ 1000$	100.00 10
att23 756 att23 756 att23 16641 att23 16641 att23 16641 att23 16641 att23 16641 att23 16621 att23 16522 att23 16522 att23 1662 att23 1552 att24 1551	754 1639 1134 1639 1137 1075 130 1845 1408 1038 1408 1038 1408 1038 1038 1038 1038 1038 1038 1038 10	$\begin{smallmatrix} 0 & 0 \\ 0 & 0 \\ 0 & 0 \\ 13 \\ 0 & 0 \\ 10 \\ 0 & 0 \\ 10 \\ 0 \\ 10 \\ 0 \\ 10 \\ 0 \\ 10 \\ 0 \\ $	o 2 11616 600 11 1122 0 2 30 1 2 31 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 9 1 1 2 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 3 3 3 3 3 3 2 3	$\begin{array}{c} 100.00\\ 90.88\\ 90.88\\ 90.88\\ 91.84\\ 91.84\\ 91.84\\ 91.84\\ 91.84\\ 91.86\\ 9$	$\begin{array}{c} 100.00\\ 99.40\\ 99.40\\ 1100.00\\ 188.40\\ 100.00\\ 99.40\\ 99.93\\ 99.93\\ 100.00\\ 99.93\\ 100.00\\ 99.93\\ 100.00\\ 99.33\\ 100.00\\ 100.00\\ 99.37\\ 99.37\\ 99.37\\ 99.37\\ 99.20\\ 100.00\\ 100.00\\ 99.37\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 99.37\\ 100.00\\ 100.00\\ 99.37\\ 100.00\\ 100.00\\ 99.46\\ 99.30\\ 100.00\\ 100.00\\ 99.46\\ 99.64\\ 100.00\\ 99.64\\ 100.00\\ 99.68\\ 99.6$
	16.39 12.44 12.44 12.44 12.07 12.65 12.55 12	$\begin{smallmatrix} 1 & 0 & 0 \\ 0 & 0 & 15 \\ 0 & 0 & 18 \\ 0 & 0 & 0 & 0 \\ 0 & 0 & 0 & 0 \\ 0 & 0 &$	2 116 6 200 19 12 2 3 3 4 4 4 3 3 3 2 2 3 8 4 4 1 1 1 2 3 3 3 2 3 8 3 2 1 1 5 1 1 1 2 3 3 8 4 4 1 1 2 3 3 2 3 8 3 2 2 1 1 2 3 3 8 4 3 2 3 2 3 3 2 3 3 2 3 3 2 3 3 3 2 3 3 3 3 2 3	09.88 09.88 01.41 01.41 01.41 01.41 01.41 01.41 01.45 0.	99.94 100.00 100.00 100.00 100.00 100.00 97.92 100.00 98.97 100.00 98.97 100.00 98.53 100.00 99.51 100.00 99.72 100.00 99.92 100.00 99.92 100.00 99.48 99.94 100.00 99.64 100.00 99.64 100.00 100.00 99.88 99.83 99.84 99.85 99.86 99.86 99.86 99.86 99.86 99.86 99.86 99.86 99.86 99.86 99.86 99.86 99.86 99.86
ek211 1247 1233 13007 ek31 752 ek31 752 ek31 752 ek31 752 ek31 1532 ek31 1532 ek31 1532 ek31 1532 ek31 1532 ek30 1532 ek30 1026 ek30 1026 ek30 1026 ek30 1026 ek30 1026 ek301 1759 ek302 1157 ek303 1026 ek304 1180 ek404 1180 ek404 1180 ek405 1031 file23 1063 file23 1063 file23 1063 file23 1064 file33 771 file33 771 file33 672 file110 572	1244 1307 1308 1308 1488 1498 1498 1498 1498 1498 1498 149	$\begin{array}{c} 0 \\ 0 \\ 1 \\ 1 \\ 1 \\ 1 \\ 1 \\ 1 \\ 1 \\ 1 \\$	116 2002 19 102 2 30 121 2 1 3 3 4 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3	$\begin{array}{c} 00.63\\ 00.63\\ 01.64\\ 01.64\\ 01.79\\ 00.78\\ 01.79\\ 00.71\\ 00$	100.00 100.00 100.00 100.00 97.92 98.97 100.00 99.93 100.00 100.00 100.00 100.00 97.32 100.00 99.87 99.87 99.87 99.87 100.00 100
algob algob algob <td>1077 1075 1075 1084 1088 1088 1088 1088 1088 1088 1098 1098 1098 1098 1098 1097 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1091 1011 1058 1061 1061 1061 1061 1061 1061 1065 1057</td> <td>$\begin{smallmatrix} 0 \\ 0 \\ 115 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\$</td> <td>6 200 1 191 0 2 30 1 1 1 3 3 4 4 3 3 2 3 3 8 4 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3</td> <td>$\begin{array}{l} 09.54\\ 09.54\\ 09.50\\ 07.50\\ 07.50\\ 09.50\\ 00$</td> <td>100.00 88.40 107.92 98.97 100.00 98.97 100.00 99.73 100.00 99.73 100.00 99.73 100.00 99.73 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.64 100.00 99.64 100.00 99.64 99.88 99.82 99.82 99.82 99.82 99.82 99.64 99.64 13.3 99.64 15.3 99.64 15.3 99.65 9</td>	1077 1075 1075 1084 1088 1088 1088 1088 1088 1088 1098 1098 1098 1098 1098 1097 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1091 1011 1058 1061 1061 1061 1061 1061 1061 1065 1057	$ \begin{smallmatrix} 0 \\ 0 \\ 115 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ $	6 200 1 191 0 2 30 1 1 1 3 3 4 4 3 3 2 3 3 8 4 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3	$\begin{array}{l} 09.54\\ 09.54\\ 09.50\\ 07.50\\ 07.50\\ 09.50\\ 00$	100.00 88.40 107.92 98.97 100.00 98.97 100.00 99.73 100.00 99.73 100.00 99.73 100.00 99.73 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.76 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.78 100.00 99.64 100.00 99.64 100.00 99.64 99.88 99.82 99.82 99.82 99.82 99.82 99.64 99.64 13.3 99.64 15.3 99.64 15.3 99.65 9
zbl31 732 zbl31 732 zbl31 8502 zbl31 8502 zbl31 1502 zbl31 1502 zbl31 1502 zbl31 1502 zbl31 1502 zbl32 1501 zbl32 1502 zbl31 1502 zbl31 1502 zbl32 1502 zbl32 1502 zbl32 1512 zbl32 1512 zbl32 1512 zbl32 1512 zbl32 1139 zbl32 1139 zbl32 1142 zbl32 1012 zbl32 1013 zbl42 1014 zbl42 1002 zbl4114 1002 zbl4114 1002 zbl4114 102 zbl4114 124 zbl4114 124 zbl4114 124	730 1 364 145 145 146 1038 1408 1038 1408 1038 1408 1038 1402 1408 1402 140 140 140 140 140 140 140 140 140 140	$\begin{array}{c} 0 \\ 118 \\ 10 \\ 1 \\ 0 \\ 10 \\ 0 \\ 0 \\ 0 \\ 11 \\ 1 \\ $	l 19 1922 0 2 301 21 21 21 21 23 3 22 3 3 22 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 2 3 2 2 3 2 3 2 2 3 2 3 2 3 2 3 2 3 2 3 2 3 2 3 2 3 3 2 2 3 3 2 3 2 3 3 2 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 3 3 2 3 3 3 3 2 3 3 3 3 3 3 2 3 3 3 3 3 2 3 3 3 3 3 2 3 3 3 3 2 3 3 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 3 3 2 3 3 3 3 2 3	$\begin{array}{c} 99 \\ 80 \\ 90 \\ 70 \\ 70 \\ 90 \\ 90 \\ 90 \\ 90 \\ 9$	$\begin{array}{c} 100.00\\ 97.92\\ 97.97\\ 99.03\\ 99.93\\ 100.00\\ 99.93\\ 100.00\\ 99.93\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 99.91\\ 99.937\\ 99.92\\ 99.937\\ 99.92\\ 99.937\\ 100.00\\ 99.937\\ 100.00\\ 100.00\\ 99.93\\ 100.00\\ 100.00\\ 99.01\\ 100.00\\ 100.00\\ 99.01\\ 100.00\\ 99.01\\ 100.00\\ 99.01\\ 100.00\\ 99.01\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.03\\ 100.00\\ 99.04\\ 100.00\\ 99.63\\ 99.63\\ 99.63\\ 100.00\\ 100.00\\ 99.63\\ 100.00\\ 1$
ktbb Stof ktb Stof ktb Stof ktb	864 1248 1248 1408 1038 1252 2006 1024 1037 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1157 1101 110 125 1205 809 1205 809 1205 809 1205 801 1006 923 755 1410 1051 801 804 804 805 806 606 607 11412 1141 1151 1152 1153 1164 1170	$\begin{array}{c} 18\\ 10\\ 0\\ 0\\ 0\\ 0\\ 0\\ 224\\ 0\\ 1\\ 2\\ 1\\ 2\\ 1\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\$	l 19 1922 0 2 301 21 21 21 21 23 3 22 3 3 22 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 3 2 2 3 2 3 2 2 3 2 3 2 2 3 2 3 2 3 2 3 2 3 2 3 2 3 2 3 2 3 3 2 2 3 3 2 3 2 3 3 2 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 3 3 2 3 3 3 3 2 3 3 3 3 3 3 2 3 3 3 3 3 2 3 3 3 3 3 2 3 3 3 3 2 3 3 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 2 3 3 3 3 2 3 3 3 3 2 3	97,50,97,97,97,97,97,97,97,97,97,97,97,97,97,	97 92 98,97 109,93 100,06 100,06 100,06 100,06 100,07
zd33 1532 zd33 1532 zd30 1234 zd30 1234 zd30 1234 zd30 1129 zd30 1139 zd30 1139 zd30 1139 zd30 1139 zd30 1139 zd30 1139 zd30 892 zd303 1112 zd303 1112 zd303 1031 L16270 1008 L16270 1008 L16270 1008 L16370 1027 L16370 1027 L16370 1028 L16370 1028 L16370 1028 L16370 1121	1265 1348 1408 1348 1408 1408 1408 1408 1408 1408 140 140 140 140 140 140 140 140 140 140	$\begin{matrix} 13\\ 0\\ 1\\ 0\\ 0\\ 0\\ 0\\ 0\\ 1\\ 2\\ 1\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\$	112 0 2 0 0 1 2 1 3 3 2 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 3 3 3 2 1 3 6 1 1 1 3 2 2 3 3 2 2 1 3 3 2 2 1 3 3 2 2 3 3 2 2 1 1 3 2 2 1 1 3 2 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 2 2 1 3 2 2 1 3 2 2 1 3 2 2 1 3 2 2 1 3 2 2 1 1 1 1	01, 79 01, 79 000, 007 07, 11 07, 11 07, 11 07, 11 07, 11 07, 11 07, 11 07, 11 09, 13 09, 13 09, 13 09, 13 09, 13 09, 23 09, 25 09, 25 09, 25 09, 25 09, 23 09, 25 09, 23 09, 25 09, 25 09, 23 09, 25 09, 23 09, 25 09, 23 09, 25 09, 23 09, 23 09, 25 09, 23 09, 25 09, 23 09, 25 09, 25 16, 6 09, 25 16, 7 10, 100 10, 10	$\begin{array}{c} 98,97\\ 100,00\\ $
action 1352 action 1352 action 1352 action 1001 action 1001 action 1001 action 1002 action 2012 action 2012 action 2012 action 2012 action 1002 action 1002 action 1002 action 1002 action 1002 action 1002 action 1159 action 1159 action 1159 action 1137 action </td <td>1348 1498 1498 1498 1498 1498 1498 1498 14</td> <td>$\begin{array}{c} 0\\ 0\\ 1\\ 0\\ 0\\ 0\\ 0\\ 24\\ 0\\ 0\\ 1\\ 2\\ 2\\ 1\\ 1\\ 1\\ 0\\ 6\\ 2\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\$</td> <td>2 3 3 1 1 2 1 3 3 4 2 9 1 3 3 2 3 3 2 3 3 2 1 3 3 2 3 2 3 2 1 1 3 3 2 2 3 2 1 1 1 3 3 2 2 3 2 1 1 1 3 3 2 2 3 2 2 1 1 1 3 3 2 2 3 2 2 3 2 1 1 1 3 3 2 2 3 2 2 3 2 2 1 1 3 3 2 2 3 3 2 2 3 2 3</td> <td>$\begin{array}{c} 100,000\\ 90,87\\ 90,87\\ 91,113\\ 80,83\\ 90,85\\ 91,80\\$</td> <td>$\begin{array}{c} 100.00\\ 99.01\\ 99.00\\ 100.00\\ 1$</td>	1348 1498 1498 1498 1498 1498 1498 1498 14	$\begin{array}{c} 0\\ 0\\ 1\\ 0\\ 0\\ 0\\ 0\\ 24\\ 0\\ 0\\ 1\\ 2\\ 2\\ 1\\ 1\\ 1\\ 0\\ 6\\ 2\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\$	2 3 3 1 1 2 1 3 3 4 2 9 1 3 3 2 3 3 2 3 3 2 1 3 3 2 3 2 3 2 1 1 3 3 2 2 3 2 1 1 1 3 3 2 2 3 2 1 1 1 3 3 2 2 3 2 2 1 1 1 3 3 2 2 3 2 2 3 2 1 1 1 3 3 2 2 3 2 2 3 2 2 1 1 3 3 2 2 3 3 2 2 3 2 3	$\begin{array}{c} 100,000\\ 90,87\\ 90,87\\ 91,113\\ 80,83\\ 90,85\\ 91,80\\$	$\begin{array}{c} 100.00\\ 99.01\\ 99.00\\ 100.00\\ 1$
Labor 1501 Abber 1502 Abber 1503 Abber 1502 Abber 1502 Abber 1507 Abber 1507 </td <td>1498 103 103 103 103 103 103 103 103 103 103</td> <td>$\begin{matrix} 1 \\ 0 \\ 0 \\ 1 \\ 1 \\ 9 \\ 0 \\ 0 \\ 1 \\ 2 \\ 1 \\ 0 \\ 1 \\ 2 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 0$</td> <td>2 300 1 21 3 849 1 3 2 3 2 3 8 3 2 1 3 3 2 3 8 3 2 1 3 3 2 3 8 3 2 1 3 3 2 3 8 4 9 1 3 3 2 3 8 4 9 1 3 3 2 3 8 4 9 1 1 1 1 1 1 1 1 1 1 1 1 1</td> <td>$\begin{array}{c} 0.98 \\ 0.91 \\ 0.91 \\ 1.11 \\ 0.91 \\ 1.11 \\ 0.91 \\ 0.$</td> <td>99 93 100,00 105,00 105,53 100,00 100,00 97,32 099,91 99,92 100,00 99,87 99,92 100,00 99,87 99,92 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 99,64 100,00 100,00 100,00 99,64 100,00 100,00 99,88 80,33 99,85 80,55 100,000 100,000</td>	1498 103 103 103 103 103 103 103 103 103 103	$\begin{matrix} 1 \\ 0 \\ 0 \\ 1 \\ 1 \\ 9 \\ 0 \\ 0 \\ 1 \\ 2 \\ 1 \\ 0 \\ 1 \\ 2 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 1 \\ 0 \\ 0$	2 300 1 21 3 849 1 3 2 3 2 3 8 3 2 1 3 3 2 3 8 3 2 1 3 3 2 3 8 3 2 1 3 3 2 3 8 4 9 1 3 3 2 3 8 4 9 1 3 3 2 3 8 4 9 1 1 1 1 1 1 1 1 1 1 1 1 1	$\begin{array}{c} 0.98 \\ 0.91 \\ 0.91 \\ 1.11 \\ 0.91 \\ 1.11 \\ 0.91 \\ 0.$	99 93 100,00 105,00 105,53 100,00 100,00 97,32 099,91 99,92 100,00 99,87 99,92 100,00 99,87 99,92 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 99,64 100,00 100,00 100,00 99,64 100,00 100,00 99,88 80,33 99,85 80,55 100,000 100,000
ck806 1040 ck807 1032 ck807 823 ck807 2012 ck808 1026 ck808 1155 ck808 1155 ck808 1155 ck809 1180 ck801 1180 ck802 4801 ck803 1026 ck804 1180 ck803 1031 ck804 1031 ck805 1031 ck804 1031 ck804 1034 ck804 1034 ck804 1034 ck804 1034 ck804 1344 ck804 1344 ck804 1344 ck804 1344 <td>1038 152 152 152 152 155 1004 101 155 1178 155 1178 155 1178 117 110 115 115 115 115 115 115 115</td> <td>$egin{array}{cccccccccccccccccccccccccccccccccccc$</td> <td>30 31 21 33 429 13 32 23 32 23 3 3 2 3 3 2 1 3 3 2 1 3 3 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 1 3 3 2 2 1 1 1 1 1 1 1 1 1 1 1 1 1</td> <td>$\begin{array}{c} 97,11\\ 99,51\\ 99,85\\ 99,85\\ 99,85\\ 99,85\\ 99,85\\ 99,85\\ 99,81\\ 99,81\\ 99,81\\ 99,81\\ 99,81\\ 99,82\\ 99,82\\ 99,82\\ 99,83\\ 99,93\\ 99$</td> <td>100.00 100.01 940.03 940.00 97.32 99.87 99.87 99.87 100.00 99.87 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.04 40.00 99.88 99.88 99.88 99.83 99.85 99.82 99.64 64.53 99.64 53.00 99.64 53.00 99.64 53.00 99.65 59.55 59.</td>	1038 152 152 152 152 155 1004 101 155 1178 155 1178 155 1178 117 110 115 115 115 115 115 115 115	$egin{array}{cccccccccccccccccccccccccccccccccccc$	30 31 21 33 429 13 32 23 32 23 3 3 2 3 3 2 1 3 3 2 1 3 3 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 3 3 2 2 1 1 3 3 2 2 1 1 1 1 1 1 1 1 1 1 1 1 1	$\begin{array}{c} 97,11\\ 99,51\\ 99,85\\ 99,85\\ 99,85\\ 99,85\\ 99,85\\ 99,85\\ 99,81\\ 99,81\\ 99,81\\ 99,81\\ 99,81\\ 99,82\\ 99,82\\ 99,82\\ 99,83\\ 99,93\\ 99$	100.00 100.01 940.03 940.00 97.32 99.87 99.87 99.87 100.00 99.87 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.03 100.00 99.04 40.00 99.88 99.88 99.88 99.83 99.85 99.82 99.64 64.53 99.64 53.00 99.64 53.00 99.64 53.00 99.65 59.55 59.
abbs 1294 abbs 1294 abbs 1294 abbs 1294 abbs 1294 abbs 1010 abbs 1195 abbs 1195 abbs 1195 abbs 1195 abbs 1197 abbs 1112 abbs 1112 abbs 1112 abbs 1111 abs	1292 2008 11157 2008 11157 21157 21157 21157 21157 21157 21157 21157 21157 2115 2115	$\begin{array}{c} 19\\ 0\\ 0\\ 24\\ 0\\ 1\\ 1\\ 2\\ 1\\ 0\\ 6\\ 2\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\$	21 3 44 29 1 3 3 2 3 3 2 3 3 2 2 3 3 2 2 3 3 2 2 1 3 3 2 1 3 3 2 1 3 4 3 1 1 1 1 1 1 2 2 1 1 1 3 4 3 1 2 2 3 2 2 1 3 5 5 5 5 5 5 5	$\begin{array}{c} 08, 33\\ 08, 33\\ 09, 33\\ 09, 34\\ 09, 36\\ 09, 74\\ 09, 36\\ 09, 74\\ 09, 36\\ 09, 74\\ 09, 36\\ 09, 36\\ 09, 36\\ 09, 36\\ 09, 30\\ 00, 30\\$	98.53 100.00 107.32 100.00 99.917 100.00 99.92 100.00 99.72 100.00 99.48 99.94 100.00 99.48 99.94 100.00 99.61 100.00 99.61 100.00 99.63 100.00 99.64 100.00 1
kalilo 2012 kalilo 2012 kalilo 1026 kalilo 1026 kalilo 1026 kalilo 1026 kalilo 1159 kalilo 1159 kalilo 1159 kalilo 1159 kalilo 1159 kalilo 113 kalilo 113 kalilo 113 kalilo 1031 kalilo 1032 kalilo 1031 kalilo 1032 kalilo 1033 kalilo 1034 kalilo 1035 kalilo 1037 kalilo 1037 kalilo 1037 kalilo 1037 kalilo 1037 kalilo 1037 kalilo 1043 kalilo 1043 kalilo 1043 kalilo 1044 kalilo	2008 1023 1157 1155 1157 1155 1157 1555 1157 1555 1157 1555 1029 105 105 105 105 105 105 105 105 105 105	$egin{array}{cccccccccccccccccccccccccccccccccccc$	3 84 29 1 3 3 2 3 8 3 2 3 8 3 2 1 3 6 1 1 1 1 5 2 1 3 28 1 1 1 1 3 28 1 1 1 3 28 2 1 3 2 2 1 3 6 1 1 3 2 2 3 6 2 1 3 6 2 3 6 2 1 3 6 2 3 6 2 3 6 2 3 6 2 3 6 2 3 6 2 3 6 2 3 6 2 3 6 2 3 6 7 2 3 6 7 2 3 6 7 2 3 6 7 2 3 6 7 2 3 6 7 3 6 7 3 6 7 3 6 7 3 6 7 3 6 7 3 6 7 3 6 7 3 6 7 3 6 7 7 1 3 6 7 7 1 1 3 6 7 7 7 1 3 7 7 7 7 7 7 7 7 7 7 7 7 7 7 7	$\begin{array}{c} 99,85\\ 99,85\\ 91,85\\ 90,95\\ 90$	$\begin{array}{c} 100, 00\\ 100, 00\\ 100, 00\\ 99, 91\\ 99, 91\\ 99, 92\\ 99, 92\\ 99, 92\\ 99, 92\\ 99, 92\\ 99, 72\\ 100, 00\\ 99, 73\\ 100, 00\\ 99, 73\\ 100, 00\\ 99, 10\\ 100, 00\\ 99, 01\\ 100, 00\\ 99, 01\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 99, 01\\ 100, 00\\ 100, 00\\ 99, 64\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 100, 00\\ 99, 64\\ 99, 63\\ 99, 64\\ 100, 00\\ 109, 64\\ 100, 00\\ 100, $
esten 3 1025 (1997) 1155 1157) 1157 1157) 1157 1157) 1157 1157) 1157 1157) 1157 1157) 1157 1157) 1157 1157 1157) 1157 115	1024 901 1157 1157 1157 919 919 915 928 929 1265 929 1265 925 1005 925 1005 925 1005 925 1005 925 1005 925 1005 925 105 105 925 105 105 925 105 105 925 105 105 105 105 105 105 105 105 105 10	$egin{array}{cccccccccccccccccccccccccccccccccccc$	84 9 1 3 2 2 3 8 3 2 1 3 3 2 1 3 4 3 5 1 1 1 5 1 1 1 3 2 3 8 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 3 2 1 1 1 3 2 1 1 1 1 3 2 1 1 1 1 3 2 1 1 1 1 3 2 1 1 1 1 3 2 1 1 1 1 3 2 1 1 1 1 1 3 2 2 1 1 1 1 1 1 1 3 2 2 1 1 1 1 1 3 2 2 1 1 1 1 1 3 2 8 1 1 1 1 1 3 2 8 2 1 1 1 1 3 2 8 3 1 1 1 3 3 3 1 1 1 3 3 3 1 1 1 3 3 3 2 5 1 1 1 3 1 1 3 2 8 3 1 1 1 3 3 3 2 5 1 1 1 3 3 3 2 5 5 5 5 5 5 5 5 5 5 5 5 5	$\begin{array}{c} 0.180\\ 0.187\\ 0.078\\ 0.$	$\begin{array}{c} 100,00\\ 70,32\\ 99,31\\ 99,91\\ 99,87\\ 99,93\\ 99,91\\ 99,87\\ 99,92\\ 99,93\\ 99,93\\ 99,93\\ 99,78\\ 90,78\\ 90,78\\ 90,78\\ 90,78\\ 90,00\\ 100,00\\ 100,00\\ 100,00\\ 100,00\\ 100,00\\ 99,03\\ 100,00\\ 99,64\\ 88\\ 89,98\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,85\\ 80,33\\ 99,64\\ 81,30\\ 99,64\\ 99,63\\ 99,65\\ 9$
setter 1 704 setter 1 704 setter 1 704 setter 1 704 setter 1 704 setter 1 704 setter 1 150 setter 1 150 sette	001 702 1157 1155 1155 1155 1155 1155 1155 115	$\begin{array}{c} 24\\ 0\\ 1\\ 2\\ 1\\ 0\\ 6\\ 2\\ 0\\ 0\\ 30\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0$	29 3 3 2 3 8 8 3 2 1 3 6 1 1 1 1 5 5 1 2 1 3 2 8 3 2 1 1 1 5 5 1 1 1 1 1 1 1 1 1 1 3 6 1 1 1 1 5 5 2 1 3 6 1 1 3 6 1 1 3 6 1 1 3 6 1 1 3 6 1 1 3 6 1 3 7 2 1 3 8 8 8 8 8 8 8 8 8 8 8 8 8 8 8 8 8 8	$\begin{array}{c} 0.5, 78\\ -0.5, 78\\ -0.5, 74\\ $	97.3.2 99.3.2 100.00 99.9.7 99.9.7 100.00 99.4 99.4 99.4 99.4 99.4 100.00 100.00 100.00 99.0 100.00 99.0 100.00 99.0 100.00 99.0 100.00 99.0 100.00 99.0 100.00 99.0 100.00 99.40 100.00 100.00 100.00 99.40 100.00
ekit1 704 wik00 1159 wik01 1159 wik01 1180 wik01 1180 wik01 1180 wik01 1180 wik01 1180 wik01 1137 wik01 1137 wik01 1137 wik01 1137 wik01 1126 wik01 1126 wik01 1126 wik01 1126 wik02 1126 wik03 1126 wik04 803 wik04 804 wik07 1126 wik011 126 wik011 1216	702 1155 1155 1155 1178 1178 1178 1178 110 110 110 105 105 105 105 105 105 105	$egin{array}{ccccc} 0\\ 1\\ 2\\ 1\\ 0\\ 6\\ 2\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\$	1 3 3 2 3 8 3 2 1 3 6 1 1 1 1 5 1 2 1 3 2 8 1 1 1 1 3 2 8 1 1 1 1 3 2 8 1 1 1 1 3 2 8 1 2 1 3 6 1 1 1 1 3 6 1 3 2 1 3 6 1 3 6 1 3 6 1 3 6 1 3 6 1 3 6 1 3 6 1 3 6 1 1 3 6 1 1 3 6 1 1 1 1	$\begin{array}{c} 0.95 \\ 0.97 \\ 0.97 \\ 0.97 \\ 0.91 \\ 0.91 \\ 0.91 \\ 0.91 \\ 0.92 \\ 0.92 \\ 0.92 \\ 0.92 \\ 0.91 \\ 0.$	$\begin{array}{c} 100.00\\ 99.07\\ 99.07\\ 99.72\\ 100.00\\ 99.78\\ 1100.00\\ 99.78\\ 100.00\\ 99.78\\ 100.00\\ 99.78\\ 100.00\\ 99.64\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 99.64\\ 100.00\\ 99.88\\ 99.83\\ 99.83\\ 99.83\\ 99.83\\ 99.83\\ 99.83\\ 99.83\\ 99.83\\ 99.64\\ 99.64\\ 99.63\\ 100.00\\ 99.64\\ 99.63\\ 100.00\\ 99.64\\ 100.00\\ 99.64\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ 99.65\\ 100.00\\ $
sk200 1159 sk210 1557 sk211 1557 sk21 1557 sk21 1557 sk23 1113 sk23 1113 sk23 1113 sk23 1113 sk23 1113 sk23 1112 sk23 1112 sk24 1003 sk25 1112 sk24 1003 sk25 1112 sk26 1112 sk26 1112 sk27 1003 sk28 1003 sk29 1112 sk29 1112 sk29 761 sk29 762 sk29	1157 1157 1157 1159 110 110 110 110 115 110 110 110 110 110	$\begin{matrix} 1\\ 2\\ 1\\ 0\\ 6\\ 2\\ 0\\ 0\\ 30\\ 0\\ 10\\ 0\\ 4\\ 0\\ 168\\ 0\\ 0\\ 0\\ 0\\ 25\\ 1\\ 121\\ 123\\ 1\\ 223\\ 1\\ 2\\ 0\\ 4\\ 83\\ 40\\ \end{matrix}$	3 3 2 3 8 3 2 1 1 3 6 1 1 1 1 5 5 1 2 1 1 5 5 1 1 1 1 3 4 3 3 2 8 3 2 2 1 1 1 3 2 8 3 2 2 1 1 3 6 1 1 1 1 5 5 1 1 1 5 5 2 1 1 1 5 5 2 1 1 1 5 5 1 1 1 5 5 1 1 1 1	$\begin{array}{c} 09, 74 \\ 09, 74 \\ 09, 62 \\ 00, 62 \\ 00, 6$	99.91 99.87 99.66 99.46 99.46 99.78 100.00 100.00 99.00 90.00 90.00 99.00 99.00 99.1 100.00 99.00 99.1 100.00 99.1 100.00 99.1 100.00 1
ed821 1557 ed840 1557 ed840 8113 ed853 8113 ed853 9117 ed853 9910 ed853 8922 ed853 8922 ed853 8922 ed853 8922 ed854 1031 fild252 1013 fild262 1031 fild262 1032 fild264 1035 fild264 1032 fild264 1032 fild264 1032 fild264 1032 fild264 1032 fild277 701 fild284 10477 fild274 1041 fild174 560 fild184 1041	1178 1796 196 196 196 197 196 196 197 196 197 196 197 196 1005 1005 1005 1005 1005 1005 1005 100	$\begin{matrix} 1 \\ 0 \\ 6 \\ 2 \\ 0 \\ 0 \\ 30 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ $	3 8 3 2 1 36 1 1 1 5 5 1 2 1 328 1 1 328 1 1 328 1 1 3 3 3 2 8 2 2 5 5 5 1	99.82 99.23 99.23 99.23 99.82 99.23 99.93 99.94 99.94 99.95 99.92 99.95 99.92 99.95 99.93 99.91 99.82 99.93 99.94 99.94 99.95 99.93 99.94 99.94 99.95 99.96 99.96 99.93 99.93 99.66 99.93 95.51 16.09.96 99.85 99.88 99.85	99 92 100.00 94.8 100.00 94.6 100.00 94.6 100.00 95.6 100.00 95.6 100.00 99.6 100.00 99.6 100.00 99.6 100.00 100.00 100.00 100.00 100.00 95.8 90.0 100.00 95.8 90.3 9 90.6 1 90.00 99.6 1 90.00 99.6 1 90.00 99.6 1 1 0 0 99.6 1 1 0 0 99.6 1 0 0 99.6 1 0 0 0 99.6 1 0 0 0 99.6 1 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0
ekki7 803 1113 113 113 113 113 113 113	7999 1111 911 913 110 913 988 857 1305 989 110 915 920 920 920 920 920 920 920 920 920 920	$egin{array}{c} 0 \\ 6 \\ 2 \\ 0 \\ 0 \\ 30 \\ 0 \\ 10 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\$	3 8 3 2 1 36 1 1 1 5 5 1 2 1 328 1 1 328 1 1 328 1 1 3 3 3 2 8 2 2 5 5 5 1	99 622 99 227 99 820 99 880 99 880 99 880 99 880 99 89 89 99 820 99 95 10 99 95 10 99 95 10 99 95 10 99 89 97 431 99 89 99 80 99 80 90 80	100.00 99.4% 99.00 90.00 90.00 90.61 90.00 99.64 99.03 99.64 99.91 100.00 99.91 100.00
ettes 3 1113 4071 2 071 4071 2 071 4072 2 071 4072 2 071 4072 0 111 4072 0 111 4072 0 111 4072 0 111 4072 0 111 4073 1 003 4074 0 103 4074	1110 1110 115 110 115 124 124 124 124 124 124 105 100 106 106 106 106 106 106 106 106 106		3 8 3 2 1 36 1 1 1 5 5 1 2 1 328 1 1 328 1 1 328 1 1 3 3 3 2 8 2 2 5 5 5 1	99 2.28 99 670 99 80 95 96 99 92 98 83 99 92 99 82 99 95 99 95 90	99.46 99.78 1100.60 96.61 100.00 99.64 100.00 99.64 100.00 99.64 100.00 100.00 100.00 100.00 100.00 100.00 99.88 99.88 99.88 99.85 99.86 99.86 99.64 99.64 99.64 99.65
cstr1 -917 cstr2 -900 cstr2 -900 </td <td>915 985 887 100 1029 1029 1005 1006 1006 1006 1006 1006 1006 1006</td> <td>$egin{array}{ccccc} 0\\ 0\\ 30\\ 0\\ 10\\ 0\\ 0\\ 4\\ 0\\ 1\\ 0\\ 168\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 25\\ 1\\ 1\\ 121\\ 23\\ 1\\ 121\\ 23\\ 1\\ 2\\ 0\\ 4\\ 83\\ 40 \end{array}$</td> <td>3 1 1 3 1 1 1 1 5 1 1 1 1 1 1 1 1 1 1 1 1 1</td> <td>.99 6.67 .99 8.0 .93 9.0 .93 9.0 .94 9.0 .95 9.0 .98 9.1 .99 9.1 .99 8.1 .99 9.8 .99 9.8 .99 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.46 .90 9.55 .16 6.00 .90 9.8</td> <td>99,78 100,00 100,00 90,00 99,00 99,00 99,00 99,00 99,00 99,00 100,00 99,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 99,88 80,33 99,85 80,55 99,55 99,64 99,65</td>	915 985 887 100 1029 1029 1005 1006 1006 1006 1006 1006 1006 1006	$egin{array}{ccccc} 0\\ 0\\ 30\\ 0\\ 10\\ 0\\ 0\\ 4\\ 0\\ 1\\ 0\\ 168\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 25\\ 1\\ 1\\ 121\\ 23\\ 1\\ 121\\ 23\\ 1\\ 2\\ 0\\ 4\\ 83\\ 40 \end{array}$	3 1 1 3 1 1 1 1 5 1 1 1 1 1 1 1 1 1 1 1 1 1	.99 6.67 .99 8.0 .93 9.0 .93 9.0 .94 9.0 .95 9.0 .98 9.1 .99 9.1 .99 8.1 .99 9.8 .99 9.8 .99 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.8 .90 9.46 .90 9.55 .16 6.00 .90 9.8	99,78 100,00 100,00 90,00 99,00 99,00 99,00 99,00 99,00 99,00 100,00 99,00 100,00 100,00 100,00 100,00 100,00 100,00 100,00 99,88 80,33 99,85 80,55 99,55 99,64 99,65
ek872 9900 ek873 859 ek873 859 ek874 850 ek874 850 ek874 850 ek874 850 ek874 850 ek874 850 ek875 850	988 857 900 1029 849 1110 1005 920 920 920 920 920 920 920 920 920 920	$egin{array}{ccccc} 0\\ 0\\ 30\\ 0\\ 10\\ 0\\ 0\\ 4\\ 0\\ 1\\ 0\\ 168\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 25\\ 1\\ 1\\ 121\\ 23\\ 1\\ 121\\ 23\\ 1\\ 2\\ 0\\ 4\\ 83\\ 40 \end{array}$	1 36 1 11 5 2 2 1 3282 1 1 1 1 1 34 3 1 1899 28 3 2 1 5 55 1	99,80 99,88 93,96 94,92 94,93 94,93 99,85 99,95 99,87 99,86 99,97 99,86 99,86 99,97 99,86 90,86 90	$\begin{array}{c} 100, 00\\ 100, 00\\ 0, 61\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\$
ektr3 859 ektr3 859 ektr3 852 ektr3 105 litter3 110 litter3 111 litter3 1008 litter3 1008 litter3 1007 litter3 1007 litter3 1007 litter3 1007 litter3 111 litter3 111 litter3 111 litter3 400 litter3 800 litter3 800 litter3 117 litter3 117 litter3 117 litter3 117 litter3 117 litter3 117 litter3 <td>857 890 1265 849 1410 1064 1064 1064 1005 1006 923 1006 923 1006 923 804 810 905 804 810 905 806 805 805 805 805 805 805 805 805</td> <td>$egin{array}{ccccc} 0\\ 30\\ 0\\ 10\\ 0\\ 4\\ 1\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 25\\ 1\\ 1\\ 121\\ 223\\ 1\\ 2\\ 0\\ 4\\ 83\\ 40\\ \end{array}$</td> <td>1 36 1 11 5 5 2 1 3282 1 1 1 1 1 34 3 1 1899 28 3 2 1 5 55 1</td> <td>99.88 95.96 98.93 98.93 99.83 99.83 99.83 99.83 99.89 67.43 99.89 99.89 99.89 99.89 99.89 99.80 99.80 99.80 99.80 99.80 99.80 99.80 99.80 99.80 72.33 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 99.80 72.53 99.80 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 90.80</td> <td>$\begin{array}{c} 100.00\\ 96.61\\ 100.00\\ 99.01\\ 109.60\\ 109.60\\ 100.00\\ 99.91\\ 100.00\\ 99.91\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 99.88\\ 99.88\\ 99.83\\ 99.83\\ 99.83\\ 99.83\\ 99.84\\ 99.84\\ 99.64\\ 99.64\\ 99.64\\ 99.64\\ 99.64\\ 99.63\\ 99.65\\ 99.$</td>	857 890 1265 849 1410 1064 1064 1064 1005 1006 923 1006 923 1006 923 804 810 905 804 810 905 806 805 805 805 805 805 805 805 805	$egin{array}{ccccc} 0\\ 30\\ 0\\ 10\\ 0\\ 4\\ 1\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 0\\ 25\\ 1\\ 1\\ 121\\ 223\\ 1\\ 2\\ 0\\ 4\\ 83\\ 40\\ \end{array}$	1 36 1 11 5 5 2 1 3282 1 1 1 1 1 34 3 1 1899 28 3 2 1 5 55 1	99.88 95.96 98.93 98.93 99.83 99.83 99.83 99.83 99.89 67.43 99.89 99.89 99.89 99.89 99.89 99.80 99.80 99.80 99.80 99.80 99.80 99.80 99.80 99.80 72.33 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 99.80 72.53 99.80 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 99.80 99.80 72.53 99.80 90.80	$\begin{array}{c} 100.00\\ 96.61\\ 100.00\\ 99.01\\ 109.60\\ 109.60\\ 100.00\\ 99.91\\ 100.00\\ 99.91\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 100.00\\ 99.88\\ 99.88\\ 99.83\\ 99.83\\ 99.83\\ 99.83\\ 99.84\\ 99.84\\ 99.64\\ 99.64\\ 99.64\\ 99.64\\ 99.64\\ 99.63\\ 99.65\\ 99.$
etess 3 522 1110222 253 1110222 253 1110223 253 111023 1112 111023 1112 11102	890 1260 1260 1260 1260 1260 1260 1260 1410 100 100 100 100 100 100 100 100 1	30 0 10 0 1 1 168 0 0 0 25 1 121 23 1 23 1 2 4 83 40 40 40 1 121 23 148 121 23 148 121 23 148 121 23 148 121 23 148 121 23 148 121 23 148 121 23 148 121 23 148 148 148 148 148 140 148 1400 1400 1400 1400 1400	36 11 11 5 1 2 1 328 1 1 1 1 1 34 3 1 1 1 1 34 3 2 8 2 8 2 1 5 1 1 1 1 1 2 1 1 2 1 1 2 1 1 2 1 1 2 1 1 2 1 1 2 1 1 2 1 1 2 1 1 2 2 1 1 2 2 1 1 2 2 1 1 2 2 2 1 1 1 2 2 2 1 1 1 1 2 2 2 1 1 1 1 1 2 2 2 1 1 1 1 1 1 1 1 1 1 1 1 1	95 596 99 023 99 023 99 055 99 055 99 889 99 869 99 90 99 90 90 90 90 90 90 90 90 90 90 90 90 90 9	96.6.1 100.00 990.64 100.00 99.64 100.00 99.01 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 99.88 99.88 99.88 99.88 99.88 99.88 99.88 99.88 99.88 99.88 99.85 99.64 100.00 99.64 100.00 99.64 100.00 99.64 100.00 99.64 100.00 99.64 100.00 99.65 100.00 99.65 100.00 99.65 100.00 99.65 100.00 99.85 99.85 99.86 100.00 99.64 100.00 99.64 100.00 99.64 100.00 99.65 100.00 100.
	1265 1029 849 1001 1001 1005 920 923 923 923 925 802 800 810 810 810 810 810 810 810 810 810	0 10 0 4 0 1 168 0 0 0 0 0 0 0 0 0 0 0 0 0	1 11 1 5 1 2 1 328 1 1 1 1 1 34 3 1 1 89 28 3 2 1 5 55 1	99 92 98 93 99 88 99 95 99 91 90 82 99 93 99 83 99 89 99 87 99 99 98 99 99 88 99 88 99 86 99 89 99 86 99 89 99 86 99 86 99 93 99 86 99 86 99 80 99 80 99 80 99 90 90 93 91 50 11 80 99 88	100.00 99.03 100.00 99.64 09.64 09.07 109.00 80.17 100.00 100.00 100.00 100.00 100.00 100.00 99.88 99.88 99.88 99.88 99.83 99.84 99.85 99.85 99.86 100.00 99.64 90.63
 Intel.46 Intel.37 Intel.38 Intel.39 Intel.38 Intel.39 Intel.38 Intel.39 Intel.38 Intel.39 Intel.38 Intel.38 Intel.39 Intel.38 Intel.38 Intel.39 Intel.38 Intel.39 Intel.38 Intel.39 Intel.30 Intel.30 Intel.30 Intel.30<	1029 849 849 849 110 1085 923 93 94 104 1085 92 93 104 81 8 81 87 0 96 66 66 65 95 112 12 141 9 41 12 141 9 41 10 41 10 11 13 14 12 14 1 14 1 14 1 14 1 14 1 1 1 1 1	10 0 4 0 1 0 168 0 0 0 0 0 225 1 121 23 1 2 0 4 83 40	11 5 1 2 1 328 1 1 1 1 1 1 1 34 3 1 189 28 3 2 1 5 793 55	98,933 99,55 99,91 99,82 99,83 99,84 99,83 99,85 99,80 90,800 90,800 90,800 90,800 90,800 90,800	99.03 100.00 99.64 100.00 90.01 90.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00 99.88 99.89 99.89 99.89 99.89 99.89 99.89 99.80 90.80 90.
 IIIICT22 851 IIIICT23 11101 IIIICT23 11102 IIIICT23 1102 IIIITT23 1102 IIIIT23 1102 IIIIT23 1102 IIIT24 IIIT25 IIIT25 IIIT25 	849 1110 1005 1200 200 1200 200 200 200 200 200 200 2	0 4 0 1 0 0 0 0 0 0 0 0 0 0 0 0 0	1 5 1 2 1 328 1 1 1 1 1 34 3 1 1 88 3 2 2 1 5 55 1	99.88 99.55 99.91 99.82 99.89 67.43 99.89 99.89 99.89 99.89 99.89 99.80 99.88 99.80 95.81 99.85 99.86 99.86 99.86 99.86 99.86 99.85 91.60 91.80	100.00 99.64 100.00 99.91 100.00 80.17 100.00 100.00 100.00 100.00 96.88 99.89 80.33 96.68 99.82 99.85 99.85 99.85 100.00 99.86 100.00 94.63 99.85 90.85 90.
111123 1112 000 000 000 000 000 000 000 000 00	1061 1085 923 923 1015 006 923 1015 002 923 923 925 923 925 923 923 924 924 924 1024 1024 1026 925 1124 925 124 925 124 925 124 124 124 124 124 124 125 124 125 125 124 125 125 125 125 125 125 125 125 125 125	$egin{array}{c} 0 \\ 1 \\ 0 \\ 168 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 0 \\ 25 \\ 1 \\ 121 \\ 23 \\ 1 \\ 2 \\ 0 \\ 4 \\ 83 \\ 40 \end{array}$	2 1 328 1 1 1 1 34 3 1 189 28 3 2 1 5 793 55 1	99.91 99.82 99.82 99.89 67.43 99.89 99.89 99.87 99.87 99.88 99.88 99.88 99.86 99.86 99.86 99.86 99.30 35.98 99.30 99.30 99.30 99.55 16.00 91.80	100.00 99.91 100.00 80.17 100.00 100.00 100.00 100.00 96.88 99.89 99.89 99.89 99.89 99.85 99.85 99.86 90.86
1016435 1087 1164539 922 116539 922 116539 922 116539 922 116738 761 117453 1047 116748 761 117453 1047 116748 761 117453 1047 116748 761 117453 1047 116748 761 117453 1047 116019 872 116111 698 116114 698 116114 698 116114 698 116113 8810 116113 8810 116113 8810 116113 1153 116120 1153 116130 880 116131 1202 116140 1102 116141 1202 116141 1202 116141 1202 116141 <	1085 920 923 923 925 925 925 925 925 925 925 925 925 925	$\begin{matrix} 1 \\ 0 \\ 168 \\ 0 \\ 0 \\ 0 \\ 0 \\ 25 \\ 1 \\ 121 \\ 23 \\ 1 \\ 2 \\ 1 \\ 2 \\ 0 \\ 4 \\ 83 \\ 40 \end{matrix}$	2 1 328 1 1 1 1 34 3 1 189 28 3 2 1 5 793 55 1	99,82 99,89 67,43 99,89 99,89 99,89 99,89 99,89 99,89 99,89 99,89 99,89 99,89 99,89 99,89 99,46 99,46 99,93 99,46 99,93 99,55 16,00 91,80 99,88	99,91 100,00 80,17 100,00 100,00 100,00 100,00 96,88 99,83 99,83 99,83 99,83 99,83 99,83 99,83 99,83 99,83 99,83 99,85 99,85 99,85 99,64 64,53 93,90 100,00 99,64
14639 922 14639 922 146795 926 146795 926 146795 926 146795 926 146795 926 146795 926 14619 927 14743 1047 14743 104	920 1006 759 759 802 804 805 806 806 806 806 806 806 807 1412 1412 1412 1412 1412 1412 1412 141	0 168 0 0 0 0 25 1 1 121 23 1 2 0 4 83 40	1 328 1 1 1 1 34 3 1 189 28 3 2 1 5 793 55 1	99.89 67.43 99.89 99.87 99.80 99.88 99.88 99.88 99.88 99.66 99.89 97.2.33 95.98 99.46 99.86 99.86 99.93 99.46 99.93 99.55 16.00 91.80 99.88	100.00 80.17 100.00 100.00 100.00 96.88 99.88 99.88 99.89 80.33 96.68 99.82 99.86 100.00 99.86 100.00 99.64 64.53 93.90 100.00
1100730 1008 1110736 025 1110736 025 1110736 025 1110736 025 1110736 025 1110736 025 1110736 025 1110736 025 1110736 025 1110736 025 11111 1111 1111 1111 1111 1111 1111	1006 923 759 1015 802 804 814 814 814 816 806 862 662 662 662 662 662 662 663 819 1119 1412 1412 1412 838 808 1119 1412 1424 1055 1535 1244 1055 1535 1244 1056 1026 867 1125 1244 1255 1255 1255 1255 1255 1255	168 0 0 0 0 25 1 1 1 23 1 2 0 4 83 40	328 1 1 1 1 3 4 3 1 1 89 28 3 2 1 5 5 793 55 1	67.43 99.89 99.87 99.89 99.89 99.89 95.81 99.66 99.89 72.33 95.98 99.46 99.89 99.46 99.86 99.93 99.55 16.00 91.80 99.88	80.17 100.00 100.00 100.00 100.00 96.88 99.88 99.88 99.80 99.80 99.86 100.00 99.66 99.86 100.00 99.64 64.53 99.64 64.53 99.64
146786 925 146786 925 146786 925 146797 97 146797 97 146797 97 146797 97 146797 97 146797 97 146797 92 14679 147 147 147 147 147 147 147 147 147 147	923 759 801 802 803 804 805 805 805 805 805 805 805 805 805 805	0 0 0 25 1 1 121 23 1 2 0 4 83 40	1 1 1 34 3 1 189 28 3 2 1 5 793 55 1	99.89 99.87 99.90 99.88 99.89 99.89 99.66 99.89 97.2.33 95.98 99.46 99.86 99.86 99.86 99.93 99.55 16.00 91.80 99.88	100.00 100.00 100.00 100.00 99.88 99.88 99.89 96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.64
116795 761 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117453 1047 117452 114 117453 1047 117452 114 117452 114 117452 114 117452 114 117452 114 117452 114 11745 114 114 114 114 114 114 114 114 114 1	759 104 85 104 85 10 87 10 87 10 87 10 87 10 87 10 87 10 11 11 11 11 11 11 11 11 11 11 11 11	0 0 0 25 1 121 23 1 2 0 4 83 40	1 1 34 3 1 189 28 3 2 1 5 793 55 1	99.87 99.90 99.88 99.89 95.81 99.66 99.89 72.33 95.98 99.46 99.89 99.46 99.93 99.46 99.93 99.55 16.00 91.80 91.80	100.00 100.00 100.00 96.88 99.88 99.89 99.82 99.82 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
117453 1047 117453 1047 117453 1047 117451 104107 823 1140107 823 1140107 823 1140107 823 1140114 650 1141 650 1141 650 1141 650 1141 650 1141 650 1141 650 1141 650 1141 650 1141 114 114 114 114 114 114 114 114 11	802 804 814 810 906 906 906 906 907 112 112 112 112 1119 944 944 944 944 944 944 94	0 25 1 1 121 23 1 2 0 4 83 40	1 1 34 3 1 189 28 3 2 1 5 793 55 1	99.88 99.89 95.81 99.66 99.89 72.33 95.98 99.46 99.86 99.93 99.55 16.00 91.80 91.80 99.88	100.00 100.00 96.88 99.89 99.89 99.89 99.89 99.80 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
kettisö 897 kettisö 897 kettisö 822 kettis 822 kettis 822 kettis 82 kettis 82 ket	894 814 814 814 814 814 814 814 814 814 84 84 84 84 84 84 84 84 84 84 84 84 84	0 25 1 121 23 1 2 0 4 83 40	1 34 3 1 189 28 3 2 1 5 793 55 1	99.89 95.81 99.66 99.89 72.33 95.98 99.46 99.86 99.93 99.55 16.00 91.80 91.80 99.88	100.00 96.88 99.89 99.89 99.89 99.80 33 99.86 100.00 99.64 64.53 93.90 100.00 99.63
kell07 823 kell07 823 kell11 67 kell11 67 kell114 658 kell114 658 kell114 658 kell114 658 kell114 1550 kell114 1414 kell124 1414 kell124 1414 kell23 045 kell23 1414 kell23 045 kell23 1414 kell23 045 kell36 810 kell36 810 kell36 810 kell36 810 kell36 1555 kell36 1555 kell36 1442 kell36 1442 kel146 1442 kell36 1442 kell36 1442 kel146 1442 k	810 870 665 666 557 1412 1112 1119 913 808 808 808 808 808 808 808 808 808 80	25 1 1 121 23 1 2 0 4 83 40	34 3 189 28 3 2 1 5 793 55 1	95.81 99.66 99.89 72.33 95.98 99.46 99.86 99.93 99.55 16.00 91.80 99.88	96.88 99.89 80.33 96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
kell 10 872 kell 11 904 kell 14 908 kell 14 143 kell 15 1414 kell 14 144 kell 15 1414 kell 15 1414 kell 15 1414 kell 15 915 kell 15 916 kell 15 1246 kell 16 1246 kell 12 1246 kell 12 1246 kell 12 1246 kell 13 1087 kell 14 127 kell 14 1088	870 906 662 957 1412 1412 1412 1413 141 941 670 838 808 1244 1661 1573 1244 1661 1573 1244 1573 1244 1575 1244 1405 1244 1244 1244 1244 1244 1244 1244 124	1 1 121 23 1 2 0 4 83 40	3 1 189 28 3 2 1 5 793 55 1	99.66 99.89 72.33 95.98 99.46 99.86 99.93 99.55 16.00 91.80 99.88	99.89 80.33 96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
belatti 908 belatti 638 belatti 638 belatti 638 belatti 638 belatti 638 belatti 1434 belatti 672 belatti 810 belatti 1246 belatti 1246 belatti 1246 belatti 1246 belatti 1246 belatti 1246 belatti 1240 belatti 1240 belatti 1240 belatti 1240 belatti 1240 belatti 1434 belatti 1434 belatti 1434 belatti 1434 <t< td=""><td>906 682 695 145 1412 1412 1412 1414 441 441 441</td><td>1 121 23 1 2 0 4 83 40</td><td>1 189 28 3 2 1 5 793 55 1</td><td>99.89 72.33 95.98 99.46 99.93 99.55 16.00 91.80 99.88</td><td>99.89 80.33 96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63</td></t<>	906 682 695 145 1412 1412 1412 1414 441 441 441	1 121 23 1 2 0 4 83 40	1 189 28 3 2 1 5 793 55 1	99.89 72.33 95.98 99.46 99.93 99.55 16.00 91.80 99.88	99.89 80.33 96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
kell12 684 kell11 684 kell11 634 kell11 1434 kell12 1414 kell13 121 kell14 121 kell13 121 kell14 121 kell12 1414 kell13 121 kell13 121 kell13 121 kell13 124 kell14 121 kell15 1246 kell16 1575 kell11 1575 kell13 1575 kell13 1246 kell13 127 kell13 127 kell13 127 kell13 127 kell14 131 kell14 131 kell14 133 kell14 136 kell14 137 kell14 137 kell14 137 kell14 137 <td>682 696 35 1412 1412 1412 944 943 808 808 808 808 808 808 808 808 808 80</td> <td>121 23 1 2 0 4 83 40</td> <td>189 28 3 2 1 5 793 55 1</td> <td>72.33 95.98 99.46 99.86 99.93 99.55 16.00 91.80 99.88</td> <td>80.33 96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63</td>	682 696 35 1412 1412 1412 944 943 808 808 808 808 808 808 808 808 808 80	121 23 1 2 0 4 83 40	189 28 3 2 1 5 793 55 1	72.33 95.98 99.46 99.86 99.93 99.55 16.00 91.80 99.88	80.33 96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
kell 14 698 kell 15 560 kell 15 560 kell 15 560 kell 15 161 kell 15 161 kell 16 112 kell 18 112 kell 18 610 kell 18 610 kell 18 610 kell 18 116 kell 18 116 kell 18 116 kell 18 126 kell 11 1022 kell 11 1028 kell 18 1038 kell 19 108 kell 10 108 kell 10 108 kell 11 1067 kell 11 1067 kell 11 1068 ke	696 557 1412 1419 640 640 640 640 640 640 640 640 640 640	23 1 2 0 4 83 40	28 3 2 1 5 793 55 1	95.98 99.46 99.86 99.93 99.55 16.00 91.80 99.88	96.68 99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
Lett 16 500 Lett 16 500 Lett 12 1434 Lett 14 1373 Lett 14 1373 Lett 14 1374 Lett 14 1374 Lett 14 1372 Lett 14 1372 Lett 14 1372 Lett 14 1374	1412 1412 1413 670 838 808 1244 1061 1573 1244 1061 1573 1244 1573 1244 1573 1244 1575 1200 1577 1200 1057 1200 1059 1440 1182 125 1200 1059 1440 1182 1395	1 2 0 4 83 40	3 2 1 5 793 55 1	99.46 99.86 99.93 99.55 16.00 91.80 99.88	99.82 99.86 100.00 99.64 64.53 93.90 100.00 99.63
Lett121 1434 Lett122 1444 Lett122 1444 Lett122 1444 Lett122 1444 Lett122 1444 Lett122 1444 Lett123 145 Lett129 672 Lett133 840 S13 Lett170 851 Lett180 1246 Lett1	1412 1412 1413 670 838 808 1244 1061 1573 1244 1061 1573 1244 1573 1244 1573 1244 1575 1200 1577 1200 1057 1200 1059 1440 1182 125 1200 1059 1440 1182 1395	2 0 4 83 40	2 1 5 793 55 1	99.86 99.93 99.55 16.00 91.80 99.88	99.86 100.00 99.64 64.53 93.90 100.00 99.63
Lieff122 1414 11 1112 112 1121 1121 1121 1121	1412 1412 441 443 4570 838 808 1064 1064 1064 1075 1075 1075 1075 1075 1075 1075 1075	0 4 83 40	1 5 793 55 1	99.93 99.55 16.00 91.80 99.88	100.00 99.64 64.53 93.90 100.00 99.63
Lield 12.4 1121 Lield 12.4 1121 Lield 12.5 042 Lield 12.5 042 Lield 12.5 042 Lield 12.5 142 Lield 13.5 142 Liel	1119 410 40 438 808 818 818 124 1573 1043 1051 1053 1043 1057 1050 1057 1000 802 957 1200 802 957 1200 802 957 1059 1059 1059 1059 1059 1059 1059 1059	83 40	793 55 1	16.00 91.80 99.88	64.53 93.90 100.00 99.63
Lett29 672 4403 440 440 440 440 440 440 440 440 44	670 838 88 811 1064 1064 1013 1214 103 1214 105 125 125 1026 1026 1026 1182 1125 1026 1182 1125 1026 1182 1125 1125 1026 1182 1125 1125 1125 1125 1125 1125 1125	40	55 1	91.80 99.88	93.90 100.00 99.63
Ledit3.3 S40 Ledit3.4 S41 Ledit3.6 S11 Ledit3.0 S12 Ledit3.0 S12 Ledit3.1 L53 Ledit3.1 L184 Ledi4.1 L196	838 800 1244 1061 1573 1244 1057 1735 1200 1029 1440 1182 507 1029 1440 1182 517 1200 1029 1440 1182 53 1040 1014 1014 1014 1014 1014 1014 101		1	99.88	100.00 99.63
ted13.6 810 ted16.0 813 ted16.0 813 ted16.0 813 ted16.0 813 ted16.0 813 ted16.0 1063 ted210 1063 ted211 1575 ted211 1075 ted214 1028 ted304 1031 ted404 1031 ted404 1032 ted504 1031 ted604 1031 ted604 1034 ted604 1034 ted604 1038 ted604 1037 ted604 1037 ted704 1045 ted33 1576 ted34 837 ted34 1366 ted34 1376 ted44 1336 ted44 1371 ted45 1337 ted46 371 ted46 3138 ted46	308 811 1244 1573 1043 1214 957 1735 1043 1214 957 1305 867 1429 1440 155 1429 1440 1155 1425 1455	0		99.88	99.63
telti6 813 telt70 897 telt70 897 telt70 126 telt70 122 telt60 142 telt60 142 telt60 1018 telt60 1014 telt70 1024 telt70 1024 telt70 1031 telt70 104 telt70 104 telt70 1054 telt7 127 telt84 1247 telt84 1247 telt84 1247 telt84 1247 telt84 1358 telt84 1358<	811 895 1244 1641 1733 1043 1214 025 1775 1026 867 1026 867 1026 1182 1125 702 1092 1014 1084 1084 1084 1084 1084 1084	3	4		
ledit70 897 ledit70 1246 ledit701 1246 ledit701 1246 ledit701 1257 ledit701 1575 ledit701 1216 ledit701 1202 ledit702 1217 ledit703 1028 ledit703 1028 ledit704 1202 ledit704 1216 ledit705 121 ledit705 122 ledit704 1242 ledit707 1184 ledit707 1184 ledit707 1184 ledit707 1184 ledit701 1018 web10 1018 web11 1016 web12 1306 web13 1307 web14 1307 web14 1307 web14 1307 web14 1308 web14 1308 web19 1398	895 1244 1063 1073 1033 1214 957 1335 1036 807 1029 1440 1185 1440 1185 1455 1056 1052 1014 1084 1084 1084 1084 1084 1084 1084	0	1	99.88	100.00
Led210 1063 Med210 11575 Med210 1575 Med210 1575 Med210 1575 Med210 1575 Med40 559 Med40 757 Med40 757 Med40 1737 Med40 1732 Med60 1031 Med60 1041 Med60 1042 Med60 1042 Med60 1042 Med60 1084 Med60 1084 Med60 1084 Med60 1127 Med60 1084 Med70 004 Med70 004 Med70 004 Med71 1062 Med31 1070 Med31 1070 Med31 1070 Med31 1070 Med31 1070 Med31 1070 Med31 1070 <	1061 1573 1043 1214 957 1735 1026 1029 1040 1029 1040 1182 1182 1182 1182 1182 1182 1182 1029 1041 1014 1044 1184 1184 1184	õ	i	99.89	100.00
Lett211 1575 Lett211 1575 Lett212 1245 Lett212 1245 Lett212 1245 Lett211 1245 Lett211 1245 Lett211 1202 Lett211 1205 Lett211 1205 <td>1573 1043 1214 957 1735 1200 1026 867 1029 1440 1185 749 1092 1092 1094 1094 1084 1194 525 895</td> <td>0</td> <td>4</td> <td>99.68 99.91</td> <td>100.00</td>	1573 1043 1214 957 1735 1200 1026 867 1029 1440 1185 749 1092 1092 1094 1094 1084 1194 525 895	0	4	99.68 99.91	100.00
led003 1045 led009 1045 led009 1216 led009 1216 led009 1236 led009 1236 led009 1028 led009 1028 led009 1028 led001 1028 led006 1442 led001 1184 led006 1442 led001 1184 led006 1444 led001 1184 led006 1404 led001 1088 led004 1088 led00 led014 1088 led00 led01	1043 1214 957 1330 1026 867 1029 1440 1182 1125 749 1092 1012 1012 1014 1044 525 895	ō	1	99.91	100.00
Leduids [216] heiden 059 heiden 059 heiden 059 heiden 1777 heiden 1777 heiden 1777 heiden 1777 heiden 1777 heiden 1787 heiden 1784 heiden 1785 heiden 17855 heiden 17855 heiden 17855 hei	1214 957 1735 1200 1026 867 1029 1440 1182 1125 749 1025 1125 749 1014 1014 1014 1014 1014 1014 1014 10	0	1	99.94	100.00
Leiduis 959 leiduin 1737 leiduin 1737 leiduin 1737 leiduin 1837 leiduin 1837 leiduin 1837 leiduin 1847 leiduin 1847 lei	957 1735 1200 1026 867 1029 1440 1182 1125 749 1092 1014 1084 1194 525 895	1	2 57	99.81 95 30	99.90 100.00
ledutyo 17337 ledutyo 2737 lettin 2022 lettin 2022 le	1735 1200 1026 867 1029 1440 1182 1125 749 1092 1014 1084 1194 525 895	0	1	95.30 99.90	100.00
Leduli 1 202 Leduli 1 202 Leduli 1 202 Leduli 1 202 Leduli 1 202 Leduli 1 1 202 Ledu	1200 1026 867 1029 1440 1182 1125 749 1092 1014 1084 1194 525 895	0	1	99.90	100.00
Letesson 10228 Letessa 869 Letessa 869 Letessa 1031 Letessa 1032 Letessa 1032 Letessa 1032 Letessa 1032 Letessa 1034 Letessa 1043 Letessa 1044 Letessa 1043 Letessa 1044 Letessa 1044 Letessa 1043 Letessa 1044 Letessa 1044 Letessa 1044 Letessa 1044 Letessa 1044 Letessa 1044 Letessa 11247 Letessa 1338 Letessa 1338 Letessa 1338 Letessa 1339 Letessa 1398 Letessa 1398 Letessa 1424 Letessa 1424 Letessa 1398 <td>1026 867 1029 1440 1182 1125 749 1092 1014 1084 1194 525 895</td> <td>0</td> <td>2</td> <td>99.83</td> <td>100.00</td>	1026 867 1029 1440 1182 1125 749 1092 1014 1084 1194 525 895	0	2	99.83	100.00
heteloidi 869 heteloidi 1031 heteloidi 1184 heteloidi 1186 heteloidi 1186 heteloidi 1186 heteloidi 1186 heteloidi 1186 heteloidi 1187 heteloidi 1187 heteloidi 1187 heteloidi 1186 heteloidi 1186 heteloidi 1186 <	867 1029 1440 1182 1125 749 1092 1014 1084 1194 525 895	0	39	96.20	100.00
hete/sh 1031 hete/sh/sh 1442 hete/sh/sh 1442 hete/sh/sh 1442 hete/sh/sh 1142 hete/sh/sh/sh/sh/sh/sh/sh/sh/sh/sh/sh/sh/sh/	1029 1440 1182 1125 749 1092 1014 1084 1194 525 895	30	84	90.33	96.32
Letter 1184 Letter 112 Letter 711 Letter 751 Letter 751 Letter 751 Letter 751 Letter 751 Letter 751 Letter 1087 Letter 1088 Letter 1087 Letter 1087 Letter 1087 Letter 1306 Letter 1366 Letter 1366 Letter 1366 Letter 1368 Letter 1368 Letter 1368 Letter 1398 Letter 749 Letter 749 Letter 749 Letter 742 Letter 742 Letter 1628	1182 1125 749 1092 1014 1084 1194 525 895	0	2	99.81	100.00
hete/sp 1127 hete/sp 1751 hete/sp 751 hete/sp 1094 hete/sp 1094 hete/sp 1094 hete/sp 1094 hete/sp 1097 hete/sp 1097 hete/sp 1171	1125 749 1092 1014 1084 1194 525 895	0	4 0	99.72	100.00 100.00
lettering 1 751 lettering 1 1094 lettering 1 1096 letteri	749 1092 1014 1084 1194 525 895	0	0 4	100.00 99.64	100.00 99.73
loga 1048 loga 1018 web3 1018 web3 1018 web3 1016 web3 1056 web3 1056 web3 527 web3 852 web3 852 web4 166 web4 164 web4 1436 web4 1437 web4 1337 web4 856 web4 856 web4 836 web4 838 web4 1398 web4 1398 web4 1398 web4 1398 web5 749 web5 749 web5 1421 17172 1628	1092 1014 1084 1194 525 895	0	4	99.64 99.87	99.73
sect0 1018 veb3 1087 veb3 527 veb3 527 veb3 527 veb3 527 veb3 527 veb3 822 veb3 822 veb3 1663 veb3 1663 veb3 1247 veb4 1366 veb4 1337 veb4 1337 veb4 1338 veb4 1339 veb4 1339 veb4 1339 veb5 749 veb5 749 veb5 1411752	1014 1084 1194 525 895	0	214	80.40	100.00
kel31 1087 wkl2 1196 wkl2 527 wkl5 822 wkl6 842 wkl6 842 wkl6 682 wkl6 1366 wkl6 1337 wkl6 1337 wkl6 1337 wkl6 931 wkl6 938 wkl6 931 wkl6 938 wkl6 938 wkl6 938 wkl6 1398 wkl6 632 wkl6 633 wkl7 856 wkl8 886 wkl7 749 wkl7 749 wkl7 1628	1084 1194 525 895	ō	3	99.70	100.00
kel22 1196 kel33 527 kel34 897 kel35 882 kel36 942 kel36 942 kel36 11663 kel40 1171 kel41 1069 kel42 1366 kel43 1247 kel44 1430 vel43 1347 kel46 971 kel47 856 kel48 838 kel48 838 kel51 749 kel51 749 kel51 749	525 895	45	385	64.52	93.96
sec14 897 sec15 882 sec16 948 sec17 682 sec18 1563 sec19 1632 sec19 1366 sec14 1136 sec14 1336 sec14 1337 sec14 1336 sec14 1337 sec14 1338 sec16 333 sec19 333 sec19 1338 sec19 1398 sec19 1394 sec19 1394 sec19 1312 sec19 1411 sec19 1411 sec19 1411	895	0	3	99.75	100.00
set35 882 set36 948 set36 948 set37 682 set48 1563 set44 1366 set44 1324 set44 1430 set44 1430 set44 1430 set45 886 set46 1337 set47 1386 set48 886 set49 1388 set51 661 set51 661 set51 1421 1732 1421 1749 1428	895	0	4	99.24	100.00
sel36 948 sel37 682 sel48 1563 sel40 1171 sel41 1069 sel42 1366 sel43 1247 sel44 1337 sel45 1337 sel46 971 sel47 856 sel48 886 sel49 1398 sel51 661 sel51 661 sel51 661 sel51 1628		0	0	100.00	100.00
sel37 682 sel38 1563 sel40 1171 sel44 1171 sel44 1337 sel45 1247 sel46 971 sel47 856 sel48 886 sel49 1398 sel51 661 sel51 661 sel51 1421 1742 1428	946	0 135	384 227	56.36 76.03	100.00 84.21
set38 1563 set40 1171 set41 1069 set42 1366 set43 1367 set44 1430 set44 1430 set45 1337 set46 971 set46 971 set47 856 set48 886 set49 1398 set51 661 set51 661 set51 1421 H7152 1628	946 679	0	511	24.74	84.21 100.00
sel40 1171 sel41 1069 sel42 1366 sel43 1247 sel44 1430 sel45 1337 sel46 971 sel47 856 sel48 856 sel49 1398 sel51 661 sel51 661 sel51 1628	1561	0	0	100.00	100.00
sch41 1069 sch42 1366 sch43 1247 sch43 1430 sch44 1430 sch46 971 sch46 971 sch46 971 sch47 856 sch48 886 sch49 1398 sch50 631 sch51 661 sch52 1411 17152 1628	1169	ő	9	99.23	100.00
sel43 1247 sel44 1430 sel45 1337 sel46 971 sel47 856 sel48 886 sel49 1398 sel45 661 sel51 749 sel52 1411 117152 1628	1067		24	97.75	100.00
sel43 1247 sel44 1430 sel45 1337 sel46 971 sel47 856 sel48 886 sel49 1398 sel45 661 sel51 749 sel52 1411 117152 1628	1363	0 2	24	98.24	99.85
xel45 1337 xel46 971 xel47 856 xel49 1398 xel50 833 xel51 661 xel51 749 xel52 1411 117152 1628	1245	0		94.94	100.00
sel46 971 sel47 856 sel48 886 sel49 1398 sel50 833 sel51 661 sel51 749 sel52 1411 147152 1628	1427		63	96.78	100.00
sel47 856 sel48 886 sel50 833 sel51 661 sel52 1411 117152 1628	1335	0	46		100.00
sel48 886 sel49 1398 sel50 833 sel51 661 sel51 749 sel52 1411 117152 1628	968	0	46 57	95.73	92.97
sel49 1398 sel50 833 sel51 661 sel51 749 sel52 1411 417152 1628	854 884	0 0 66	46 57 96	90.09	100.00
sel50 833 sel51 661 sel51 749 sel52 1411 417152 1628	884 1396	0 0 66 0	46 57 96 98	90.09 88.52	
sel51 661 sel51 749 sel52 1411 417152 1628	1396 831	0 0 66 0 0	46 57 96 98 88	90.09 88.52 90.05	100.00
sel51 749 sel52 1411 417152 1628	659	0 0 66 0 0 0	46 57 96 98 88 4	90.09 88.52 90.05 99.71	100.00 100.00
sel52 1411 117152 1628	747	0 0 66 0 0 0 0	46 57 96 98 88 4 4	90.09 88.52 90.05 99.71 99.52	100.00 100.00 100.00
117152 1628	1409	0 0 66 0 0 0 0	46 57 96 98 88 4 4	90.09 88.52 90.05 99.71 99.52 95.14	100.00 100.00 100.00 100.00
	1626	0 0 66 0 0 0 0 0 0	46 57 96 98 88 4 4 32 29	90.09 88.52 90.05 99.71 99.52 95.14 96.12	100.00 100.00 100.00 100.00 100.00
114046 1260	1258	0 0 66 0 0 0 0 0 0 0 0	46 57 96 98 88 4 4 4 32 29 1	90.09 88.52 90.05 99.71 99.52 95.14 96.12 99.93	100.00 100.00 100.00 100.00 100.00 100.00
114157 1081	1079	0 66 0 0 0 0 0 0 0 0 0 0 0	46 57 96 98 88 4 4 32 29 1 0 0	90.09 88.52 90.05 99.71 99.52 95.14 96.12 99.93 100.00 100.00	100.00 100.00 100.00 100.00 100.00 100.00 100.00 100.00
114172 663 115814 1036	661 1034	0 0 66 0 0 0 0 0 0 0 0 0 0	46 57 96 98 88 4 4 32 29 1 0	90.09 88.52 90.05 99.71 99.52 95.14 96.12 99.93 100.00	100.00 100.00 100.00 100.00 100.00 100.00 100.00

Publication	Method	# Beats	Annotation (File Name)	\mathbf{SE}	$+\mathbf{P}$
This work	Blocks of interest	111,201	Automatic (.pu)	95.0	98.59
This work	Blocks of interest	3542	Manual (.q1c)	98.90	98.77
Martinez et al. [38]	Wavelet	3542	Manual (.q1c)	99.77	97.79
Laguna <i>et al.</i> [39]	Low-pass-differentiator	3542	Manual (.q1c)	99.0	97.74
Vila $et al.$ [40]	Modelling	3542	Manual (.q1c)	96.2	N/R

Table 4: T waves detection performance comparison on the QT database [24]. (N/R: not reported).

4 Limitations of Study and Future Work

The preliminary results are promising, especially after testing the algorithm on the QT database; however, testing the algorithm on a larger sample size is necessary to generalize the findings. In addition, a more focused study is needed to investigate atrial fibrillation, atrial flutter, paroxysmal supraventricular tachycardia, junctional tachycardia, and multifocal atrial tachycardia as the morphology of T waves may differ.

The presented method assumes that the R peaks are correctly detected. There is a linear correlation between the detection of T waves and the detected R peaks. If the R peaks are misclassified, this method will fail as it depends on the position of R peaks. However, this study provides a positive proof of concept for detecting T waves in arrhythmic ECG beats.

We created a rule-based system based on the three RR interval window proposed in [32]. Perhaps, it is important to investigate the development of a rule-based system based on more than three RR intervals to examine if it will improve the overall T wave detection accuracy. Moreover, the proposed method is based on several thresholds, which are calibrated using clinical knowledge. Clinical knowledge is considered the gold standard in our work; thus, the obtained thresholds were considered optimal. However, an open question for future work is to optimize all thresholds to improve accuracy.

One of the next steps regarding the results of this study is to detect arrhythmic ECG beats using the RT or ST interval as a main feature. In addition, the detection of P waves based on the accurate detection of T wave peaks needs to be examined. Moreover, perhaps, an optimization over the clinical parameters after splitting the databases into a training set and test set may improve the detection rate of the T waves. There is also a need to investigate the T-waves with different morphology, e.g., biphasic T-waves.

Technically, exploring the event-related moving average methodology for detecting events in ECG signals is promising in terms of computational complexity and efficiency. This can be further improved by investigating other bandpass filters with different orders and also by developing fast-moving average (or median) techniques for real-time analysis and mobile phone applications.

5 Conclusions

There is a limitation when evaluating T wave detection algorithms as datasets with annotated T waves are lacking. Consequently, comparing existing algorithms becomes even more difficult. Therefore, annotation of T waves is discussed and provided. Using these approaches, it is possible to support diagnostic analysis, delivering important information for (differential) diagnosis to medical experts. Using mobile technologies with automatic analysis software driven by medical expert knowledge, it will be furthermore possible to provide screening and monitoring solutions in places where medical expertise is scarce, such as remote rural areas and developing countries.

The use of two moving averages is simple and computationally efficient for mobile electronic health tools, such as cell phones and telemedicine technologies. The assessment of the T detector has been reliably carried out over the existing standard databases (QT and MIT-BIH), which contain different beat types and morphologies found in ECG signals. The developed algorithm was evaluated on all ECG recordings in the MIT-BIH database, 48 self-annotated records containing a total of 109,985 heart beats. It achieved a sensitivity of 95% and a positive predictivity of 98.59% over the MIT-BIH ECG signals, which contain low SNR, baseline wander, paced beats, and various arrhythmias. Interestingly, the proposed algorithm succeeded in scoring the highest overall accuracy of 98.84% over the manually annotated QT database (3542 heart beats) when compared to the other three algorithms (cf. Table 4). Moreover, the algorithm scored a high overall accuracy of 96.7% over the automaticlly annotated QT database (111,201 heart beats). Overall, simplicity and efficiency are required in developing T wave detection algorithms for processing long-term recordings and large databases as well as for expanding our telemedicine capabilities in the near future.

6 Acknowledgments

Mohamed Elgendi gratefully acknowledges the Australian government and Charles Darwin University whose generous scholarships towards his PhD set the foundation for this research. He appreciates the support of Friso De Boer and acknowledges Gari Clifford for helpful discussions.

7 Author Contributions

ME designed the experiment. ME, BE, and DA performed the statistical analysis. ME, BE, and DA conceived of the study and drafted the

manuscript. The authors approved the final manuscript.

8 Conflicts of Interest

The authors declare no conflict of interest.

References

- Alwan A (2011) Global status report on noncommunicable diseases
 2010. World Health Organization .
- Dilaveris PE, Gialafos EJ, Sideris SK, Theopistou AM, Andrikopoulos GK, et al. (1998) Simple electrocardiographic markers for the prediction of paroxysmal idiopathic atrial fibrillation. American Heart Journal 135: 733–738.
- Tsipouras MG, Fotiadis DI, Sideris D (2002) Arrhythmia classification using the RR-interval duration signal. In: Proc. IEEE Computers in Cardiology. pp. 485–488.
- Tran T, McNames J, Aboy M, Goldstein B (2004) Prediction of paroxysmal atrial fibrillation by analysis of atrial premature complexes. IEEE Trans Biomedical Engineering 51: 561–569.
- de Chazal P, O'Dwyer M, Reilly RB (2004) Automatic classification of heartbeats using ECG morphology and heartbeat interval features. IEEE Trans Biomedical Engineering 51: 1196–1206.
- Krasteva VT, Jekova II, Christov II (2006) Automatic detection of premature atrial contractions in the electrocardiogram. In: *Elec*trotechniques Electronics E & E. volume 9–10, pp. 49–55.

- Scholz UJ, Bianchi AM, Cerutti S, Kubicki S (1997) Vegetative background of sleep: spectral analysis of the heart rate variability. Physiol Behavior 62: 1037–1043.
- Trinder J, Kleiman J, Carrington M, Smith S, Breen S, et al. (2001) Autonomic activity during human sleep as a function of time and sleep stage. J Sleep Res 10: 253–264.
- Zapanta L, Poon C, White D, Marcus C, Katz E (2004) Heart rate chaos in obstructive sleep apnea in children. In: Proc. 26th Annual International Conference of the IEEE Engineering in Medicine and Biology Society (IEMBS '04). pp. 3889–3892.
- Shouldice R, O'Brien L, O'Brien C, de Chazal P, Gozal D, et al. (2004) Detection of obstructive sleep apnea in pediatric subjects using surface lead electrocardiogram features. Sleep 27: 784–792.
- Elgendi M, Eskofier B, Dokos S, Abbott D (2014) Revisiting QRS detection methodologies for portable, wearable, battery-operated, and wireless ECG systems. PLOSONE 9: e84018.
- Goutas A, Ferdi Y, Herbeuval JP, Boudraa M, Boucheham B (2005)
 Digital fractional order differentiation-based algorithm for P and Twaves detection and delineation. ITBM-RBM Elsevier 26: 127-132.
- Wan X, Xu D (2010) An ECG T waves detection scheme based on the compensatory criterion. In: 2010 3rd International Conference on Biomedical Engineering and Informatics (BMEI). volume 2, pp. 730-734. doi:10.1109/BMEI.2010.5640074.
- Vazquez-Seisdedos C, Neto J, Maranon Reyes E, Klautau A, Limao de Oliveira R (2011) New approach for t-wave end detection on electro-

cardiogram: Performance in noisy conditions. BioMedical Engineering OnLine 10: 77, DOI: 10.1186/1475-925X-10-77.

- Nair A, Marziliano P (2014) P and t wave detectionon multichannel ecg using fri. In: Engineering in Medicine and Biology Society (EMBC), 2014 36th Annual International Conference of the IEEE. pp. 2269-2273. doi:10.1109/EMBC.2014.6944072.
- Clifford GD, Azuaje F, McSharry P (2006) Advanced Methods and Tools for ECG Data Analysis. Norwood, MA, USA: Artech House, Inc.
- Hayden G, Brady W, Perron A, Somers M, Mattu A (2002) Electrocardiographic T-wave inversion: Differential diagnosis in the chest pain patient. The American Journal of Emergency Medicine 20: 252– 262.
- Michael S, Brady W, Perron A, Mattu A (2002) The prominent T wave: Electrocardiographic differential diagnosis. The American Journal of Emergency Medicine 20: 243–251.
- Smith D, Nowacki D, Li JJ (2010) Ecg t-wave monitor for potential early detection and diagnosis of cardiac arrhythmias. Cardiovascular Engineering 10: 201-206.
- Moody GB, Mark RG (2001) The impact of the MIT-BIH arrhythmia database. IEEE Engineering in Medicine and Biology Magazine 20: 45–50.
- Goldberger AL, Amaral LAN, Glass L, Hausdorff JM, Ivanov PC, et al. (2000) PhysioBank, PhysioToolkit, and PhysioNet: Components

of a new research resource for complex physiologic signals. Circulation 101: e215–e220.

- 22. Willems JL, Arnaud P, Bemmel JHV, Bourdillon PJ, Degani R, et al. (1987) A reference data base for multilead electrocardiographic computer measurement programs. Journal of the American College of Cardiology 10: 1313–1321.
- 23. Elgendi М (2015)Annotation of Р and Т waves in MIT-BIH Arrhythmia Database. Available online: http://www.elgendi.net/databases.htm (accessed on 20 July 2015).
- Laguna P, Mark RG, Goldberg A, Moody GB (1997) A database for evaluation of algorithms for measurement of QT and other waveform intervals in the ECG. In: Proc. IEEE Computers in Cardiology, Sept 7–10, 1997, Lund, Sweden, pp. 673–676.
- Elgendi M (2013) Fast QRS detection with an optimized knowledgebased method: Evaluation on 11 standard ECG databases. PLOSONE 8: e73557.
- Elgendi M, Jonkman M, De Boer F (2010) Frequency bands effects on QRS detection. In: the 3rd International Conference on Bio-inspired Systems and Signal Processing (BIOSIGNALS2010), January 20–23, 2010, Valencia, Spain, pp. 428–431.
- Elgendi M, Norton I, Brearley M, Abbott D, Schuurmans D (2013) Systolic peak detection in acceleration photoplethysmograms measured from emergency responders in tropical conditions. PLOSONE 8: e76585.

- Elgendi M, Norton I, Brearley M, Abbott D, Schuurmans D (2014)
 Detection of a and b waves in the acceleration photoplethysmogram.
 BioMedical Engineering OnLine 13: 139.
- Elgendi M (2014) Detection of c, d, and e waves in the acceleration photoplethysmogram. Computer Methods and Programs in Biomedicine 117: 125–136.
- Thakor NV, Webster JG, Tompkins WJ (1983) Optimal QRS detector. Medical and Biological Engineering 21: 343-50.
- Sahambi JS, Tandon SN, Bhatt RKP (1997) Using wavelet transforms for ECG characterization. An on-line digital signal processing system. IEEE Engineering in Medicine and Biology Magazine 16: 77–83.
- Roskamm H, Csapo G (1982) Disorders of cardiac function. New York: Dekker, 1st edition.
- 33. Elgendi M, Jonkman M, De Boer F (2009) Recognition of T waves in ECG signals. In: Proc. IEEE 35th Annual Northeast Bioengineering Conference, Boston, MA, USA. pp. 1-2, DOI:10.1109/NEBC.2009.4967754.
- 34. Trahanias P (1993) An approach to QRS complex detection using mathematical morphology. IEEE Transactions on Biomedical Engineering 40: 201–205.
- 35. Laguna P, Thakor N, Caminal P, Jané R, Yoon HR, et al. (1990) New algorithm for QT interval analysis in 24-hour Holter ECG: performance and applications. Medical and Biological Engineering and Computing 28: 67–73.

- 36. Schimpf R, Wolpert C, Bianchi F, Giustetto C, Gaita F, et al. (2003) Congenital short QT syndrome and implantable cardioverter defibrillator treatment: inherent risk for inappropriate shock delivery. Journal of Cardiovascular Electrophysiology 14: 1273–1277.
- 37. Braunwald E, Zipes D, Libby P, Bonow R (2004) Braunwald's Heart Disease: A Textbook of Cardiovascular Medicine. Philadelphia: Saunders, 7th edition.
- 38. Martínez JP, Almeida R, Olmos S, Rocha AP, Laguna P (2004) A wavelet-based ECG delineator: evaluation on standard databases. IEEE Transactions on Biomedical Engineering 51: 570–581.
- 39. Laguna P, Jané R, Caminal P (1994) Automatic detection of wave boundaries in multilead ECG signals: validation with the CSE database. Computers and Biomedical Research 27: 45–60.
- 40. Vila J, Gang Y, Presedo J, Fernandez-Delgado M, Barro S, et al. (2000) A new approach for TU complex characterization. IEEE Transactions on Biomedical Engineering 47: 764–772.